

Instructions for use

RITE

Alta2
Nera2
Ria2
Alta
Nera
Ria



oticon
life-changing technology

Model overview

This booklet is valid for Oticon RITE in the following hearing aid families and models:

Inium Sense

- | | | |
|------------------------------------|------------------------------------|-----------------------------------|
| <input type="checkbox"/> Alta2 Pro | <input type="checkbox"/> Nera2 Pro | <input type="checkbox"/> Ria2 Pro |
| <input type="checkbox"/> Alta2 | <input type="checkbox"/> Nera2 | <input type="checkbox"/> Ria2 |

Inium

- | | | |
|-----------------------------------|-----------------------------------|----------------------------------|
| <input type="checkbox"/> Alta Pro | <input type="checkbox"/> Nera Pro | <input type="checkbox"/> Ria Pro |
| <input type="checkbox"/> Alta | <input type="checkbox"/> Nera | <input type="checkbox"/> Ria |

The following speakers are available for the above models:

- Speaker 60
- Speaker 85
- Speaker 100
- Power flex mould speaker 100
- Power flex mould speaker 105 (only Inium Sense models)

Introduction to this booklet

This booklet guides you in how to use and maintain your new hearing aid. Please read the booklet carefully including the **Warning section**. This will help you to achieve the full benefit of your new hearing aid.

Your hearing care professional has adjusted the hearing aid to meet your needs. If you have additional questions, please contact your hearing care professional.

For your convenience this booklet contains a navigation bar to help you navigate easily through the different sections.

| **About** | [Start up](#) | [Handling](#) | [Options](#) | [Warnings](#) | [More Info](#) |

Intended use

The hearing aid is intended to amplify and transmit sound to the ear and thereby compensate for impaired hearing within mild to severe-to-profound hearing loss. The hearing aid is intended for use by adults and children older than 36 months.

IMPORTANT NOTICE

The hearing aid amplification is uniquely adjusted and optimised to your personal hearing capabilities during the hearing aid fitting performed by your hearing care professional.

Table of content

About

Identify your hearing aid speaker and earpiece	8
Your hearing aid	10
Identify left and right hearing aid	12
Battery	13

Start up

Turn the hearing aid ON and OFF	14
When to replace a battery	15
Replace the battery	16

Handling

Caring for your hearing aid	18
MultiTool for handling batteries and cleaning	19
Put on the hearing aid	20
Standard earpieces: dome or Grip Tip	22
Replace dome or Grip Tip	23
ProWax miniFit filter	24
Replace ProWax miniFit filter	25


	Customised earpieces:	
	Micro mould, LiteTip or Power mould	26
	Clean Micro mould, LiteTip or Power mould	27
	Replace ProWax filter in customised earpieces	28
Options		
	Optional features and accessories	30
	Mute the hearing aid	31
	Change programs	32
	Change volume	34
	Tamper-resistant battery drawer	36
	Wireless accessories	38
	Other options	40
Warnings		
	Warnings	42
More info		
	Troubleshooting guide	48
	Water resistant	50
	Conditions of use	51
	Technical information	52
	Warranty	63

Identify your hearing aid speaker and earpiece

This will make it easier for you to navigate through this booklet.

Speakers

RITE hearing aids come with three different speakers:.

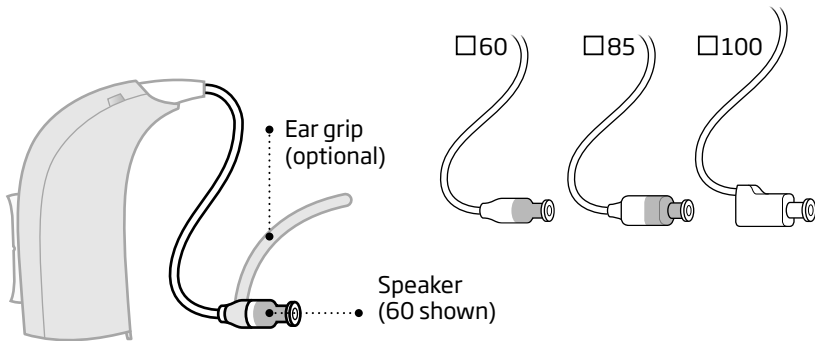

Earpieces

The speaker uses one of the following earpieces:



Dome (open dome shown)*



Grip Tip*



Micro mould or LiteTip

Power flex mould speakers


100
 105




* Please see details for replacing the dome or Grip Tip in section: Replace dome or GripTip

Your hearing aid

What it is

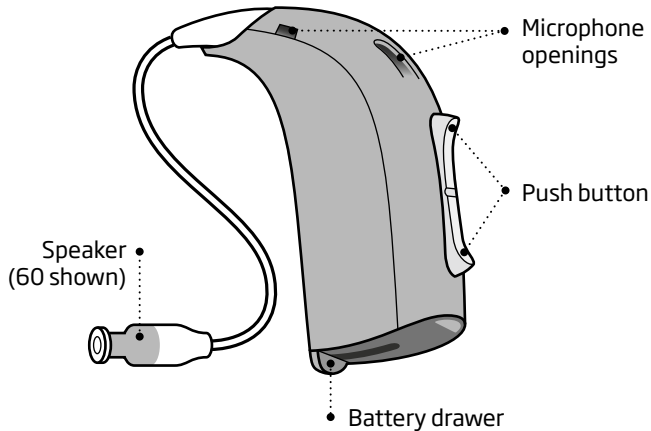

What it does

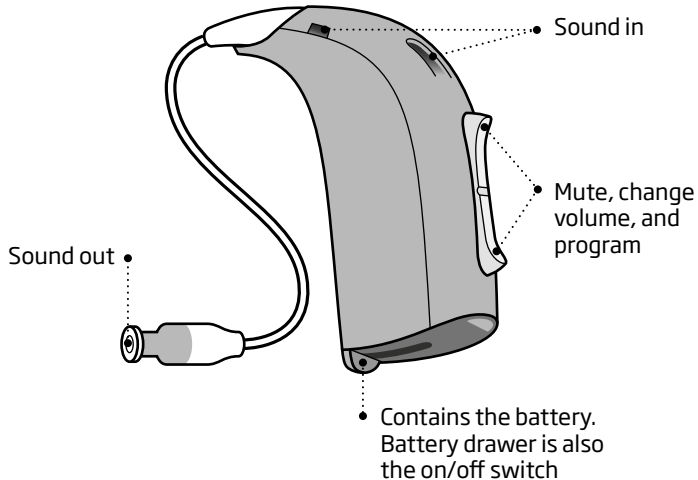


Identify left and right hearing aid

It is important to distinguish between the left and the right hearing aid, as they might be programmed differently.

You can find left/right colour indicators in the battery drawer and on 60 and 85 speakers as shown. Indicators can also be found on 100 speakers and some earpieces.


A **RED** indicator marks the
RIGHT hearing aid.

A **BLUE** indicator marks the
LEFT hearing aid.

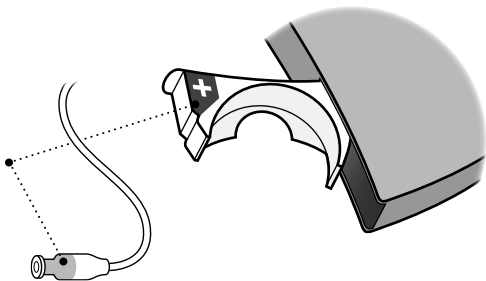

Battery

Your hearing aid is a miniature electronic device that runs on special batteries. To activate the hearing aid, you must insert a new battery in the battery drawer. See how in chapter “Replace the battery”.

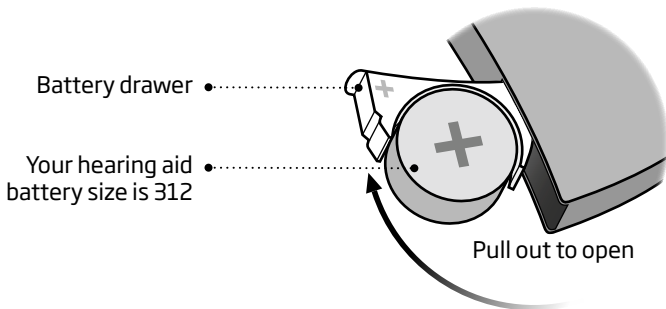


Turn the hearing aid ON and OFF

The battery drawer is also used to switch the hearing aid on and off. To preserve the battery, make sure your hearing aid is switched off when you are not wearing it.


Turn ON

Close the battery drawer with the battery in place.

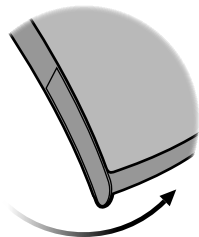

Turn OFF

Open the battery drawer.

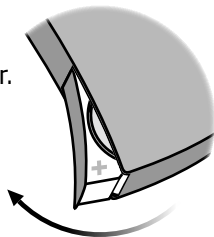

When to replace a battery

When it is time to replace the battery you will hear two beeps repeated in moderate intervals until the battery runs out.



Two beeps

= The battery is running low



Four beeps


= The battery has run out

Battery maintenance tip

To make sure the hearing aid is always working, bring spare batteries with you, or replace the battery before you leave home.


Replace the battery

1. Remove

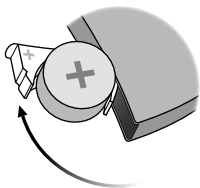

Fully open the battery drawer. Remove the battery.


2. Uncover

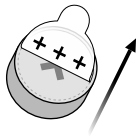

Remove the sticky label from the + side of the new battery.


Tip:
Wait 2 minutes so that the battery can draw air, to make sure to function optimally.

3. Insert

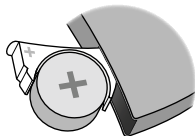


Insert the new battery into the battery drawer. Make sure the + side faces up.

4. Close

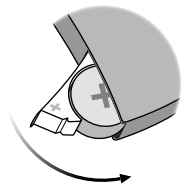

Close the battery drawer. The hearing aid will play a jingle through the earpiece. Hold the earpiece close to your ear to hear the jingle.

Tip

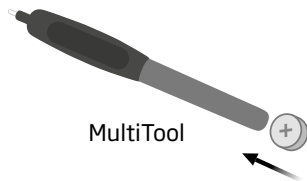

The MultiTool can be used for battery change. Use the magnetic end to remove and insert batteries.


The MultiTool is provided by your hearing care professional.

Caring for your hearing aid

When handling your hearing aid, hold it over a soft surface to avoid damage if you drop it.

Clean the microphone openings

Carefully brush away debris from the openings. Gently brush the surface. Make sure the brush is clean and that it is not pressed into the openings.

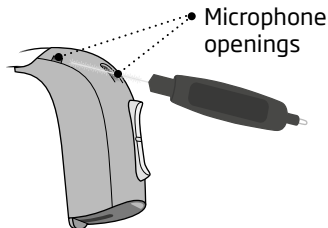


IMPORTANT NOTICE

Use a soft, dry cloth to clean the hearing aid. It must never be washed or immersed in water or other liquids.

MultiTool for handling batteries and cleaning

The MultiTool contains a magnet which makes it easier to replace the battery in the hearing aid. It also contains a brush and a wire loop for cleaning and removing ear wax. If you need a new MultiTool please contact your hearing care professional.

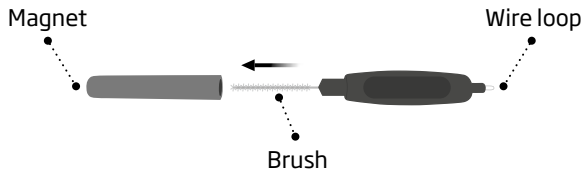


IMPORTANT NOTICE

The MultiTool has a built-in magnet. Keep the MultiTool 30 cm away from credit cards and other magnetically sensitive devices.

Put on the hearing aid

The speaker brings the sound into your ear. The speaker should always be used with an earpiece attached. Use only parts designed for your hearing aid.

If the speaker has an ear grip, place it in the ear so it follows the contour of the ear (see step 3).

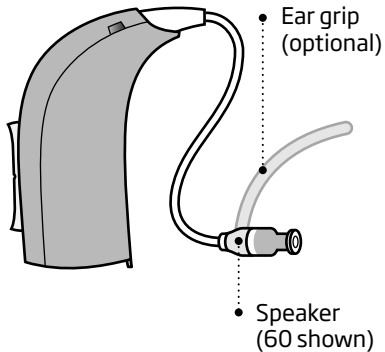


Step 1



Place the hearing aid behind your ear.

Step 2

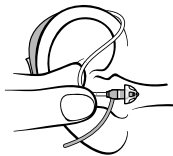

Hold the bend of the speaker wire between your thumb and index finger. The earpiece should point towards the ear canal opening.

Step 3



Gently push the earpiece into your ear canal until the speaker wire sits close against the side of your head.

Standard earpieces: dome or Grip Tip

Both domes and Grip Tip are made of soft rubber material. There are 4 different types of domes. Check your dome type and size below.



Open dome




Bass dome, Single vent



Power dome



Bass dome, Double vent

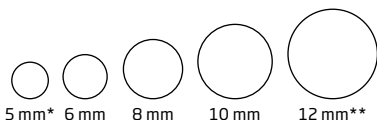

GripTip

Available in small and large, left and right with and without vent.

* only as open dome for speaker 60

** not for open dome


IMPORTANT NOTICE

If the dome is not on the speaker when removed from the ear, the dome might still be in the ear canal. Consult your hearing care professional for further instructions.

Replace dome or Grip Tip


The earpiece (dome or Grip Tip) should not be cleaned. If the earpiece is filled with wax, replace it with a new one. Grip Tip should be replaced at least once a month.

Step 1

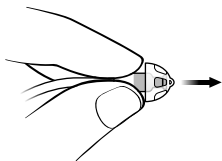


Hold on to the speaker and then pull off the earpiece.

Step 2

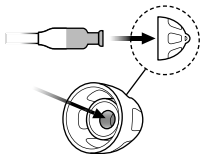

Insert the speaker exactly into the middle of the earpiece to obtain a firm attachment.

Step 3

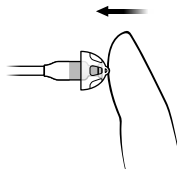

Push firmly to ensure that the earpiece is fastened securely.

ProWax miniFit filter

The speaker has a white wax filter attached to the end where the earpiece is attached. The wax filter avoid wax and debris from damaging the speaker. If the filter becomes clogged, please replace the wax filter or contact your hearing care professional.


Remove the earpiece from the speaker before replacing the wax filter.

IMPORTANT NOTICE

Always use the same type of wax filter as was originally supplied with the hearing aid. If you are in any doubt about the use or replacement of wax filters, contact your hearing care professional.


Replace ProWax miniFit filter

1. Tool

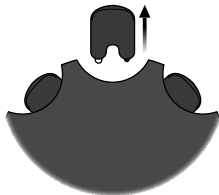


Remove the tool from the shell. The tool has two pins, one empty for removal and one with the new wax filter.

2. Remove

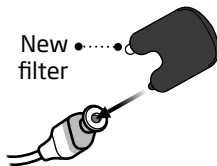

Push the empty pin into the wax filter in the speaker and pull it out.

3. Insert

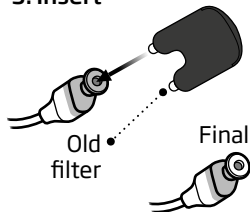

Insert the new wax filter pin using the other pin, remove the tool and throw it out.

Note:


If you use a mould it requires that the hearing care professional replaces the wax filter in the speaker.

Customised earpieces: Micro mould, LiteTip or Power mould

There are three different types of moulds: Micro mould, LiteTip, and Power mould. The earpieces are customised for your ear.



• Micro mould

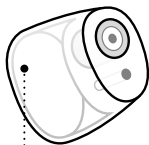

• LiteTip




• Power mould

Clean Micro mould, LiteTip or Power mould

The earpiece should be cleaned regularly.


The earpiece contains a wax filter which should be replaced when clogged or when the hearing aid does not sound normal.

Always use the wax filter delivered by your hearing care professional.

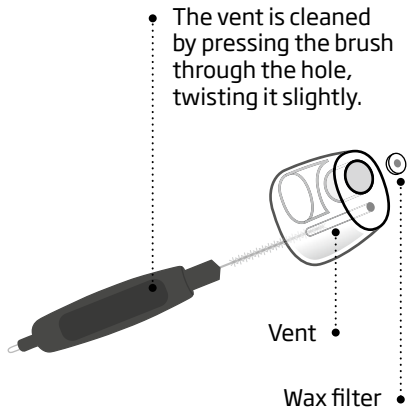


Replace ProWax filter in customised earpieces

1. Tool

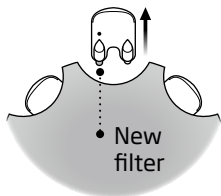


Remove the tool from the shell. The tool has two pins, one empty for removal and one with the new wax filter.

2. Remove

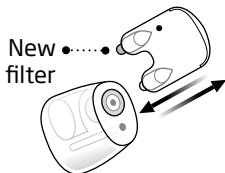

Push the empty pin into the wax filter in the micro mould and pull it out.

3. Insert

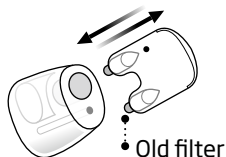

Insert the new wax filter using the other pin, remove the tool and throw it out.

Optional features and accessories


The features and accessories described on the following pages are optional. Please contact your hearing care professional to find out how your hearing aid is programmed.

If you experience difficult listening situations, a special program may be helpful. These are programmed by your hearing care professional.

Write down hearing situations where you may need help.

Mute the hearing aid

Use the mute function if you need to silence the hearing aid while wearing it.

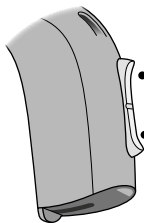


Apply a very long press to either end of the push button. To reactivate the hearing aid, push the button briefly.

IMPORTANT NOTICE

Do not use the mute function as an off switch, as the hearing aid still draws current from the battery in this mode.

□ **Change programs**



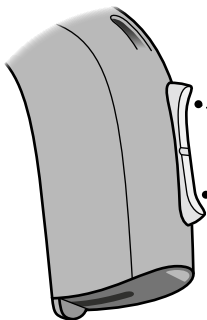

Your hearing aid can have up to 4 different programs. These are programmed by your hearing care professional.



Press up or down between programs

Note that you can change continuously between programs - both up and down in the program order. If for example you want to go from program 1 to 4, you can press the down button once instead of pressing the up button 3 times.


To be filled out by the hearing care professional

Program	Sound you will hear when activated	When to use
1	 "1 beep"	
2	 "2 beeps"	
3	 "3 beeps"	
4	 "4 beeps"	


Program change:	<input type="checkbox"/> LEFT	<input type="checkbox"/> RIGHT	<input type="checkbox"/> Short press	<input type="checkbox"/> Long press
-----------------	-------------------------------	--------------------------------	--------------------------------------	-------------------------------------

Change volume

The push button allows you to adjust the volume.
You will hear a click when you turn up or down the volume.

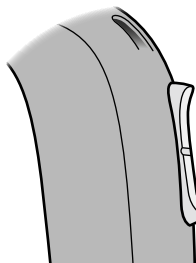

- • Press button to turn up volume.
- • Press button to turn down volume.

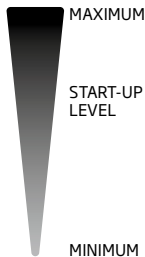

You will hear
2 beeps at
start-up level


To be filled out by the hearing care professional

Volume change	<input type="checkbox"/> LEFT	<input type="checkbox"/> RIGHT	<input type="checkbox"/> Short press
----------------------	-------------------------------	--------------------------------	--------------------------------------


Quick reset

If you wish to return to the standard settings of the hearing aid programmed by your hearing care professional, simply open and then close the battery drawer.

Open

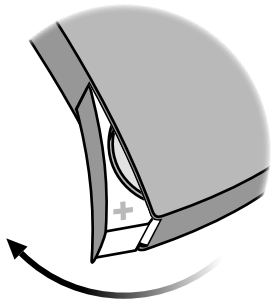

Close

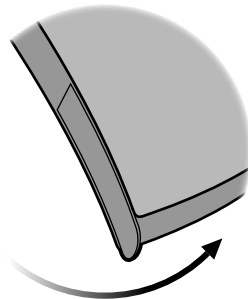

Tamper-resistant battery drawer


To keep the battery out of reach of infants, small children and people with learning difficulties, a tamper-resistant battery drawer should be used. Use a small screwdriver to open the drawer.

IMPORTANT NOTICE

Avoid using excessive force to open the battery drawer. Do not force the battery drawer beyond its fully opened position. Make sure to insert the battery correctly.


Please contact your hearing care professional if you suspect the tamper-resistant effectiveness and locking ability is damaged.

Unlock the battery drawer

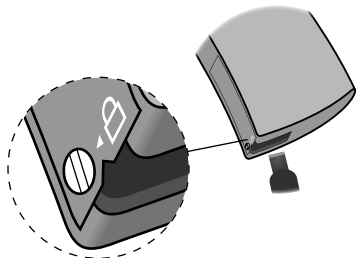

Use the screwdriver to turn the screw to unlock position (vertical). Insert the screwdriver into the slit in the bottom of the hearing aid and push the door open.

Lock the battery drawer

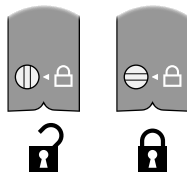

Make sure the drawer is completely closed. Turn the screw to the locked position (horizontal).

Wireless accessories

As an enhancement to your wireless hearing aid a broad range of wireless accessories are available.

ConnectLine

ConnectLine is a family of products that allows you to receive audio signals from TVs, phones, music players, PCs or an external microphone wirelessly through your hearing aid.

Remote Control

Remote Control offers the ability to change program, adjust volume or mute your hearing aid.

□ **FM**

Your hearing aid supports the option of attaching an FM system. The FM system allows you to receive speech or audio signals directly to your hearing aid from a dedicated transmitter.


For detailed information on wireless accessories contact your hearing care professional or visit **www.oticon.com.au**

Other options



Telecoil

helps you hear better when using a telephone with a built in loop or when you are in buildings with teleloop systems such as theatres, churches or lecture rooms. This symbol or a similar sign is shown wherever a teleloop has been installed.

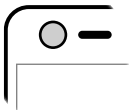

Autophone

can automatically activate a phone program in the hearing aid, if your telephone has a dedicated magnet. The magnet needs to be placed at your telephone next to the sound outlet.

□ **Direct Audio Input (DAI)**

allows your hearing aid to receive signals directly from external sources such as TV, radio, music players, etc. The DAI adaptor is mounted to your hearing aid and via a cable connected to an external sound source.

For detailed information on these options contact your hearing care professional or visit **www.oticon.com.au**

IMPORTANT NOTICE

When the DAI is connected to an audio source plugged into a wall outlet, the audio source must comply with IEC-60065, IEC-60601 or equivalent safety standards.

Warnings

You should familiarise yourself fully with the following general warnings before using your hearing aid for your personal safety and to ensure correct use.

Please note that a hearing aid will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions. Furthermore, note that in most cases, infrequent use of a hearing aid does not permit a user to attain full benefit from it.

Consult your hearing care professional if you experience unexpected operations or events with your hearing aid.

Usage of hearing aids

- Hearing aids should be used only as directed and adjusted by your hearing care professional. Misuse can result in sudden and permanent hearing loss.
- Never allow others to wear your hearing aid as incorrect usage could cause permanent damage to their hearing.

Choking hazards & risk of swallowing batteries

- Hearing aids, their parts, and batteries should be kept out of reach of children and anyone who might swallow these items, or otherwise cause injury to themselves.

- Batteries have occasionally been mistaken for pills. Therefore, check your medicine carefully before swallowing any pills.
- Most hearing aids can be supplied with a tamper-resistant battery drawer upon request. This is strongly recommended for infants, small children and people with learning difficulties.
- Children younger than 36 months must always use a tamper-resistant battery drawer. Please talk to your hearing care professional about the availability of this option.

If a battery or hearing aid is swallowed, see a doctor immediately.

Battery use

- Always use batteries recommended by your hearing care professional. Batteries of low quality may leak and cause bodily harm.
- Never attempt to recharge your batteries and never dispose of batteries by burning them. There is a risk that the batteries will explode.

Warnings

Dysfunction

- Be aware of the possibility that your hearing aid may stop working without notice. Keep this in mind when you depend on warning sounds (e.g. when you are in traffic). The hearing aids may stop functioning, for instance if the batteries have expired or if the tubing is blocked by moisture or ear wax.

Active implants

- Caution must be taken with active implants. In general, follow the guidelines recommended by manufacturers of implantable defibrillators and pacemakers regarding use with mobile phones and magnets.

- The Autophone magnet and MultiTool (which has a built-in magnet) should be kept more than 30 cm away from the implant, e.g. do not carry it in a breast pocket.
- If you have an active brain implant, please contact the manufacturer of your implantable device for information about the risk of disturbance.

X-ray, CT, MR, PET scanning and electrotherapy

- Remove your hearing aid before X-ray, CT/MR/PET scanning electrotherapy, surgery etc. as your hearing aid may be damaged when exposed to strong fields.

Heat and chemicals

- The hearing aid must never be exposed to extreme heat e.g. left inside a parked car in the sun.
- The hearing aid must not be dried in microwave ovens or other ovens.
- The chemicals in cosmetics, hairspray, perfume, aftershave lotion, suntan lotion and insect repellent can damage the hearing aid. Always remove your hearing aid before applying such products and allow time to dry before use.

Power instrument

- Special care should be exercised in selecting, fitting and using a hearing aid where maximum sound pressure capability exceeds 132 dB SPL (IEC 711), as there may be risk of impairing the remaining hearing of the hearing aid user.

For information on whether your hearing aid is a power instrument, see the back of this leaflet.

Warnings

Possible side effects

- Hearing aids, moulds or domes may cause an accelerated accumulation of ear wax.
- The otherwise non-allergenic materials used in hearing aids may in rare cases cause a skin irritation or other side effects.

Please seek consultation with a physician if these conditions occur.

Interference

- The hearing aid has been thoroughly tested for interference, according to the most stringent international standards. However, interference between the hearing aid and other devices (e.g. some mobile telephones, citizens band systems and shop alarm systems and other devices) may occur. If this occurs, increase the distance between the hearing aid and the interfering device.



WARNING

New and used button/coin batteries are hazardous and should be kept away from children. Button/coin batteries can cause severe or fatal injuries in 2 hours or less if swallowed or placed inside any part of the body. Seek medical attention immediately if this is suspected.

If it is suspected a button/coin battery has been swallowed or otherwise placed inside any part of the body, you should contact the

Australian Poisons Information Centre
on 13 11 26

New Zealand National Poisons Centre
on 0800 764 766
immediately for 24/7 fast, expert advice.

Troubleshooting guide

Symptom	Possible causes
No sound	Worn-out battery
	Clogged tube or earpieces (dome, GripTip, or mould)
Intermittent or reduced sound	Clogged sound outlet
	Moisture
	Worn-out battery
Squealing noise	Hearing aid earpiece inserted incorrectly
	Ear wax accumulated in ear canal

If none of the above solutions work, consult your hearing care professional for assistance.

Solutions

Replace the battery

Clean mould

Consider replacing wax filter, dome or Grip Tip

Clean mould or replace wax filter, dome or Grip Tip

Wipe battery with a dry cloth

Replace the battery

Re-insert hearing aid

Have ear canal examined by your doctor

Water resistant

Your hearing instrument is water resistant which means it is designed to be worn in all daily life situations. Therefore you do not have to worry about sweat or getting wet in a rain shower. Should your hearing instrument come in contact with water and stop working, please follow these guidelines:

1. Gently wipe off any water
2. Open the battery drawer and remove the battery and gently wipe off any water in the battery drawer
3. Let the hearing instrument dry with the battery drawer left open for approximately 30 minutes
4. Insert a new battery

IMPORTANT NOTICE

Do not wear your hearing instrument while showering or participating in water activities. Do not immerse your hearing instrument in water or other liquids.

Conditions of use

Operating conditions	Temperature: +1°C to +40°C Relative humidity: 5% to 93%, non-condensing
Storage and transport conditions	Temperature and humidity should not exceed the following limits for extended periods during transportation and storage: Temperature: -25°C to +60°C Relative humidity: 5% to 93%, non-condensing

Technical information

The hearing aid contains a radio transmitter using short range magnetic induction technology working at 3.84 MHz. The magnetic field strength of the transmitter is $< -42 \text{ dB}\mu\text{A/m @ 10m}$.

The emission power from the radio system is well below international emission limits for human exposure. For comparison, the radiation of the hearing aid is lower than unintended electromagnetic radiation from, for example, halogen lamps, computer monitors, dishwashers, etc. The hearing aid complies with international standards concerning Electromagnetic Compatibility.

Due to the limited space available on the hearing aids, all relevant approval markings are found in this document.

USA and Canada

The hearing aid contains a module with:

FCC ID: U28FU2BTERIT
IC: 1350B-FU2BTERIT

The device complies with Part 15 of the FCC Rules and with Industry Canada's licence-exempt RSSs.

Operation is subject to the following two conditions:

1. This device may not cause harmful interference.
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

More information: www.oticon.com.au

This Class B digital apparatus complies with Canadian ICES-003.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined

by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the manufacturer or an experienced radio/TV technician for help.

The manufacturer declares that this hearing aid is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU.

Declaration of Conformity is available from the manufacturer.



Manufactured by:
Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark
www.oticon.global

CE 0543



Waste from electronic equipment must be handled according to local regulations.



IP58

Australian Sponsor

Audmet Australia Pty Ltd
Suite 4, Level 4, Building B
11 Talavera Road
North Ryde NSW 2113

Description of symbols used in this booklet



Warnings

Text marked with a warning symbol must be read before using the device.



Manufacturer

The device is produced by the manufacturer whose name and address are stated next to the symbol. Indicates the medical device manufacturer, as defined in EU Directives 90/385/EEC, 93/42/EEC and 98/79/EC.



CE mark

The device complies with Medical Device Directive 93/42/EEC. The four digits number indicates the identification of the notified body.



Electronic waste (WEEE)

Recycle hearing aids, accessories or batteries according to local regulations. Hearing aid users can also return electronic waste to their hearing care professional for disposal. Electronic equipment covered by Directive 2012/19/EU on waste and electrical equipment (WEEE).



Regulatory Compliance Mark (RCM)

The device complies with electrical safety, EMC and radio spectrum requirements for devices supplied to the Australian or New Zealand market.

IP58

IP code

Indicates the class of protections against harmful ingress of water and particulate matter according to EN 60529:1991/A1:2002. IP5X indicates dust protection. IPX8 indicates the protection against the effects of continuous immersion in water.

Description of symbols used on the regulatory packaging label



Keep dry

Indicates a medical device that needs to be protected from moisture.



Caution symbol

Consult instructions for use for warnings and cautions.



Catalogue number

Indicates the manufacturer's catalogue number so that the medical device can be identified.



Serial number

Indicates the manufacturer's serial number so that a specific medical device can be identified.



Button Battery Warning

The device is intended to be used with button batteries. Refer to enclosed warnings before using the device.

International warranty

Your device is covered by an international limited warranty issued by the manufacturer for a period of 12 months from the date of delivery to you. This limited warranty covers manufacturing and material defects in the device itself, but not in accessories such as batteries, tubing, speakers, earpieces and filters, etc. Problems arising from improper/incorrect handling or care, excessive use, accidents, repairs made by an unauthorised party, exposure to corrosive conditions, physical changes in your ear, damage due to foreign objects entering the device, or incorrect adjustments are NOT covered by the limited warranty and may void it. Under this warranty, the

manufacturer will choose, at its sole discretion, whether to repair the device, or replace it with an equivalent model. The above warranty does not affect any legal rights that you might have under applicable national legislation governing sale of consumer goods. Your hearing care professional may have issued a warranty that goes beyond the clauses of this limited warranty. Please consult him/her for further information.

If you need service

Take your device to your hearing care professional, who may be able to sort out minor problems and make adjustments immediately. Your hearing care professional can assist you to obtain warranty service from the manufacturer. Your hearing care professional may charge a fee for their services.

Australia only

The rights given under this warranty are in addition to other rights and remedies under a law in relation to the good or services to which the warranty relates to. Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Australian customers should contact Oticon Australia at Suite 4, Level 4, Building B, 11 Talavera Road, North Ryde NSW 2113, telephone 1800 633 738 if you have any questions regarding this warranty.

Warranty

Certificate

Name of owner: _____

Hearing care professional: _____

Hearing care professional's address: _____

Hearing care professional's phone: _____

Purchase date: _____

Model left: _____ Serial no.: _____

Model right: _____ Serial no.: _____

Power instrument Yes No

Settings overview for your hearing aid

Left

Right

Yes No

Volume control

Yes No

Yes No

Program shift

Yes No

Yes No

Mute

Yes No

Volume control indicators

On Off

Beeps at min /max volume

On Off

On Off

Clicks when changing volume

On Off

On Off

Beeps at preferred volume

On Off

Battery indicators

On Off

Low battery warning

On Off



oticon
life-changing **technology**