A close-up photograph of a Bernafon hearing aid and its wax guard resting on a piano keyboard. The hearing aid is a small, beige, behind-the-ear device with a thin tube leading to a clear wax guard. The piano keys are white and black, and the lighting is soft, creating a warm and elegant atmosphere.

TURN A LITTLE TUNE INTO PURE HARMONY

For a First Class
Hearing Experience

JUNA 9 | 7
PICO RITE

bernafon[®]
Your hearing · Our passion





Bernafon is delighted to introduce the new Juna Pico RITE (receiver-in-the-ear) hearing aid in the premium performance categories.

Juna Pico RITE is incredibly small and stylish, yet powerful and comfortable to wear. With the addition of impressive Audio Efficiency™ features and the new miniFit speaker and dome system, you can offer Juna to even more clients. Plus, Juna Pico RITE is compatible with all Bernafon wireless accessories, including the SoundGate Mic and SoundGate App, in connection with the SoundGate 3.

**Make more out of every single sound.
With Juna Pico RITE. Your first choice.**

Turn a Small-Size Instrument into a Big-Impact Solution



1 Battery Door

- Easy to open and close
- Integrated stopper for turning off the hearing instrument, preventing the battery from falling out

2 Push Button

- Large, easy-to-find push button
- Program and volume changes possible via push button
- Mute possible via push button

3 Microphones

- Integrated dual microphone system supporting optimal directional performance

4 Durability


- Hydrophobic coating of key components repels liquids
- IP57 dust and moisture protection

5 Wireless Compatibility

- Connects to all Bernafon wireless accessories, such as the SoundGate 3 Bluetooth® interface

6 Programming Port

- Discreet placement of the programming slot
- Easy to insert FlexConnect programming strip

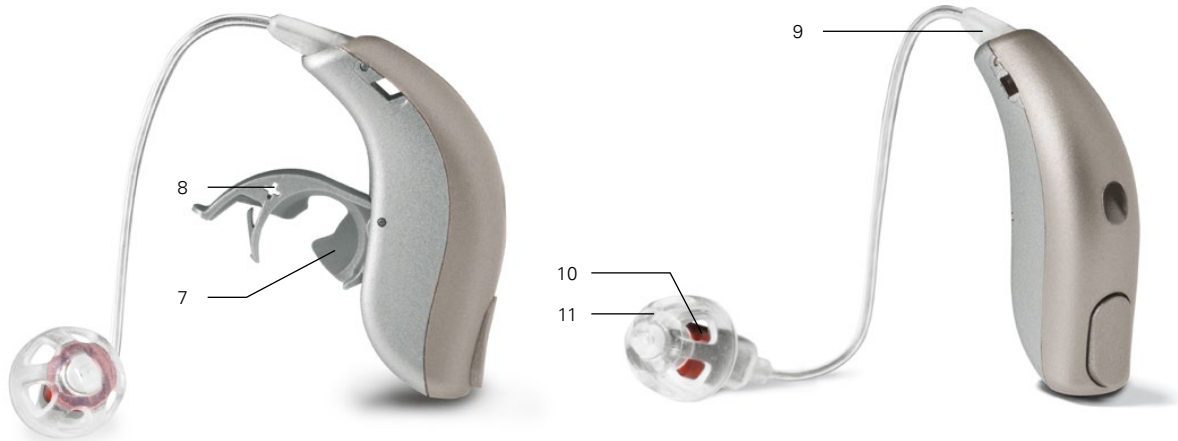

7 Battery Compartment

- Extra small battery compartment designed for the long-lasting 312 battery

8 Left/Right Indicator

- Visual indicator helps to distinguish between the left and right hearing instrument

9 miniFit Speaker Unit

- Easy-to-mount, click-on speaker unit system
- Three different speaker types for mild to severe hearing losses
- Four different speaker lengths

10 miniFit Receiver

- Visual indicator helps to distinguish between the left and right speaker unit
- Exchangeable wax filter protects the receiver from wax and moisture

11 miniFit Earpieces

- Wide range of different miniFit domes in multiple sizes
- Wide range of different custom molds

Turn Fitting Options into Real Benefits

Juna Pico RITE is incredibly small, yet powerful enough to fit hearing losses of up to 105 dB HL. It comes with the new miniFit speaker and dome system providing major benefits in fitting flexibility, performance and appearance, as well as reliability.

Fitting Flexibility

Four different speaker types are available for mild to profound hearing losses. Additionally, various solutions for instant and custom fittings are available, increasing the range of options for you and your clients.

Performance and Appearance

The 60- and 85-Speaker come in a curved housing designed to match the curve of the ear canal to support best sound quality and physical fit. The 100-Speaker has a transparent housing for discreetness. All receiver types are as small as possible, allowing for even smaller custom molds.

Reliability

The miniFit domes are made out of silicon and provide a high retention force on the receiver and in the ear. This new material helps to strengthen the durability of the domes. Integrated wax filters in the custom molds prolong the longevity of the earpieces too.

	60-SPEAKER	85-SPEAKER	100-SPEAKER	105-SPEAKER
INSTANT FIT				
Open Dome (6, 8, 10 mm)	●	●	—	—
Bass Dome, Double Vent (6, 8, 10, 12 mm)	●	●	●	—
Bass Dome, Single Vent (6, 8, 10, 12 mm)	●	●	●	—
Power Dome (6, 8, 10, 12 mm)	●	●	●	—
CUSTOM FIT				
Micro Mold	●	●	—	—
Lite Tip	●	●	—	—
Power Mold	—	—	●	●

Turn Challenging Situations into Success Stories



Meticulously developed high-end features are packed into the tiny Pico RITE hearing instrument. All emanate from Bernafon's proprietary Audio Efficiency™ technology. This core system improves your clients' speech understanding and enhances their listening comfort. Individual lifestyle needs and preferences can also be easily accommodated.

Audio Efficiency™

The SoundGate 3 is our Bluetooth® interface that allows clients to connect to multiple audio sources. The RC-N remote control combines high tactility with small size. Also introduce your clients' to the easy-to-use SoundGate App. Volume and program changes on the Pico RITE can simply be done via an iPhone® ¹⁾.

Find out more about Bernafon's latest technology in the new Audio Efficiency™ Flash application available on our website.

**RC-N
Remote
Control**



**TV Adapter 2
Phone Adapter 2**



SoundGate Mic



SoundGate 3



SoundGate App



Apple, the Apple logo, and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc.

¹⁾ The SoundGate App is for use with an iPhone (SoundGate App supports iPhone 6 Plus, iPhone 6, iPhone 5s, iPhone 5c, iPhone 5, iPhone 4s, and iPhone 4). The SoundGate App requires Bernafon's SoundGate 3 (or newer versions) and wireless hearing aids.

Turn State-of-the-Art Technology into a Comprehensive Care Package

	JUNA 9 PICO RITE	JUNA 7 PICO RITE
SIGNAL PROCESSING		
ChannelFree™	●	●
Speech Cue Priority™	●	●
Frequency Composition™	●	●
Frequency Bandwidth	10 kHz	10 kHz
LISTENING COMFORT		
Adaptive Noise Reduction Plus (ANR Plus)	5 ctr	4 ctr
Transient Noise Reduction	●	●
Adaptive Feedback Canceller Plus	●	●
Wind Noise Monitor	●	●
Environment Optimizer	1/4	1/1
Advanced Soft Noise Management	●	●
Reverb Reduction	●	-
i-VC	●	●
BINAURAL COORDINATION		
Volume Control, Program Change	●	●
Environment Classification	●	●
Non-Telephone Ear Attenuation (Auto-T)	●	●
ENTERTAINMENT		
Live Music Program	●	●
Cinema Program	●	●
DIRECTIONALITY		
Adaptive Directionality	●	●
True Directionality™	●	-
CONVENIENCE FEATURES		
VC Clicks	●	●
Mute Via Push Button	●	●
Configurable Start-Up Delay	●	●
CUSTOMIZATION		
Program Options/Memories	14/4	13/4
Adaptivity Control	●	-
Data Logging & Data Learning	●	●
VC Learning	●	-
Smart VC	●	-
Language Specific Targets	●	●
REMfit™	●	●
Client Interactive	●	●
Comfort in Airplane Program	●	-

Manufacturer:
 Bernafon AG
 Morgenstrasse 131
 3018 Bern
 www.bernafon.com

Local Manufacturer and Distributor:
 Bernafon Canada
 500 Trillium Drive, Unit 15
 Kitchener, Ontario
 Canada N2R 1A7
 www.bernafon.ca

www.bernafon.com

bernafon[®]
 Your hearing • Our passion