

Obersee, Rapperswil-Jona, Switzerland



CLEAR NATURAL SOUND



VÉRITÉ

True performance matters

www.bernafon.com

bernafon[®]
Your hearing • Our passion

True performance matters

We are surrounded by an auditory scene with multiple textures and subtle details. This soundscape provides us with a vital connection to other people and to our surroundings. Hearing loss disrupts this connection and prevents us from fully interacting with others through the world of sound.

To truly compensate for hearing loss, a hearing system must provide an accurate reproduction of the auditory world. “Vérité” is the French word for truth. Bernafon’s Vérité stays true to the sound source. Adapting automatically to changes in the sound environment, Vérité ensures that you won’t miss a thing.

Vérité’s organic design reflects the natural world. Wireless communication between both hearing instruments creates a realistic sound experience. Vérité also learns and adapts to your preferences for volume in different listening situations. Vérité delivers when true performance matters.



PERFORMANCE



DESIGN



CONNECTIVITY



High precision performance

True sound quality

For hearing instrument users, sound quality is of the utmost importance. To truly capture the sound scene, sound quality must be high for both speech and environmental sounds.

Precision is the key to Vérité's true sound experience. The power of Bernafon's proprietary ChannelFree™ digital signal processing is combined with outstanding multi-environment performance. Advanced adaptive functionality allows Vérité to adjust precisely to changes in the sound environment. This winning partnership delivers an authentic listening experience.

True signal processing

When you watch an old movie, you immediately notice that the colors seem faded, the transitions are less smooth, and the overall quality is much less sharp and clear. Just as films have made amazing progress, hearing instrument signal processing has evolved to unprecedented levels.

Vérité's signal processing stays true to the sound source. Slower-acting systems can miss the short consonants and important transitions between the different elements of speech. Vérité analyzes incoming sounds thousands of times per second, ensuring that every sound is amplified precisely. The result is clear, natural sound and an effortless listening experience.

BERNAFON VÉRITÉ

PRECISION



PERFORMANCE





BERNAFON VÉRITÉ
BALANCE

A matter of balance

Binaural coordination

When your hearing is normal, both ears act together as a system. Slight differences between the level and frequency content of sounds reaching the two ears provide spatial awareness and the ability to perceive the depth of the sound environment.

Until recently, hearing instruments have acted independently, creating an artificial listening situation. Vérité uses wireless communication to connect the instruments, restoring a natural link between the two ears.

Easy volume and program changes

Despite the sophistication of adaptive hearing systems, it is sometimes necessary to make changes to the volume or the listening program. If you wear two hearing aids, this can be tiresome when it is necessary to make changes to both instruments independently.

With Vérité's wireless coordination, changes made on one instrument are automatically transferred to the other. This makes operation of your instruments more convenient, accurate and efficient.

Smarter anti-whistling functionality

Feedback cancellation does a good job of preventing whistling in hearing instruments but sometimes certain musical tones can be suppressed by the feedback system. Vérité's intelligent feedback cancellation system is excellent for listening to music. Using binaural communication, feedback detection is coordinated between the instruments. This means that if a feedback-like signal is detected in both instruments at the same time, the system will not regard it as feedback, allowing you to enjoy music to the fullest.

PERFORMANCE





PERFORMANCE



Adapts to any environment

Your hearing needs change constantly during the course of a day. Whether you are indoors or out, in quiet or noise, conversing with one person or many, your hearing instruments must perform well in any situation.

Vérité has an advanced multi-environment program that adapts automatically to most listening situations. Sophisticated detectors constantly monitor the sound environment and make fine adjustments whenever needed to provide smooth transitions between listening situations.

Vérité instantly detects:

- Conversation in quiet
- Conversation in noise
- Noise
- Quiet
- Wind

Vérité's detection systems are coordinated to ensure that both hearing instruments work together to bring you the best hearing in every situation.

Conversation is the focus

Vérité can focus on either conversation or environmental sounds. Whenever speech is detected, the hearing system automatically changes the settings in both instruments to ensure that conversation has the highest priority. In the absence of speech, Vérité focuses on surrounding sounds.

Dynamic contrast

Communication is especially difficult when there is loud background noise. For this situation, Vérité has a special solution.

When Vérité detects speech in situations with loud background noise, this causes the hearing instruments to automatically shift to the Dynamic Contrast Mode.

In this mode, Vérité ensures that loud speech sounds are amplified clearly. At the same time, sophisticated features like Adaptive Noise Reduction and Adaptive Directionality work together to maintain comfort in background noise.

Once the loud background noise has decreased, the instruments fade back from the Dynamic Contrast Mode to normal operation.

BERNAFON VÉRITÉ

ADAPTIVITY

Panoramic sound

Vérité makes smooth transitions between speech, noise, quiet and all combinations between to accurately represent the surrounding sound scene. The sound picture is enhanced by communication between the two instruments, keeping sounds in the correct balance and rebuilding the natural partnership between the two ears.

Thanks to advanced technology, you can experience effortless listening that is customized for the situations where you most need hearing assistance. Vérité brings clear, natural sound to every situation.



PERFORMANCE



Dedicated listening programs

There are some special listening situations that require slightly different operation of the features in the hearing instruments. Vérité offers dedicated programs you can activate with the push button for these specific listening situations.

These include individual programs for television, live music, car travel and more. Your hearing care professional will assist you in choosing the programs that will be most useful for your individual lifestyle.



BERNAFON VÉRITÉ
FIDELITY



BERNAFON VÉRITÉ
STYLE

Style and design

Perfect profile

Vérité's rounded, organic form reflects harmony with nature. The user-friendly shape is easy to handle and fits very comfortably behind the ear.

Ergonomic accessories

Optional Vérité accessories include the choice of the RC-P remote control or the SoundGate communication device. Both units provide ergonomic design for easy and discreet operation of the hearing system.



SoundGate




RC-P remote control

Colors

Vérité offers a choice of many attractive colors, taken from the natural world. You are sure to find one that suits your individual preferences.



DESIGN

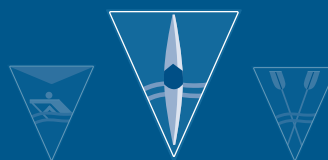



BERNAFON VÉRITÉ
WIRELESS

Wireless connectivity



SoundGate is a wireless communication link that provides a connection between your hearing instruments and audio sources, allowing the hearing instruments to work like wireless headphones.

Vérité can receive wireless sound via the SoundGate from cellular phones and other audio sources, such as TV equipped with a Bluetooth® transmitter.

As an alternative, television or other sound equipment may be directly connected by an audio cable to the SoundGate.

Whatever the source, the SoundGate then transmits the signal to your Vérité instruments wirelessly, allowing you to enjoy your favorite entertainment devices.

SoundGate advantages

- Hands free cell phone connection
- Interference-free operation
- Standby mode ready to receive calls
- Binaural reception of calls
- Remote operation of the hearing instrument
- Enhanced TV enjoyment

CONNECTIVITY



At Bernafon, we are passionate about developing quality hearing systems that enable our customers to enjoy authentic hearing experiences. Since 1946, Bernafon has been at the forefront of hearing system technology. Thanks to Swiss Engineering and precision technology, Bernafon has pioneered many market firsts including the first digitally programmable hearing instrument and our proprietary ChannelFree™ signal processing. Today, Bernafon representatives and employees in more than 70 countries fulfil our vision to help people with impaired hearing to communicate without limitation.

Bernafon AG

Morgenstrasse 131
3018 Bern
Switzerland
Phone +41 31 998 15 15
Fax +41 31 998 15 90
www.bernafon.com

Bernafon, LLC

200 Cottontail Lane, Bldg. B
Somerset, NJ 08873
USA
Phone +1 888-941-4203
Fax +1 732-560-4877
www.bernafon-us.com

Bernafon Canada Ltd.

500 Trillium Drive Unit 15
Kitchener, Ontario N2R 1A7
Canada
Toll free 1 800 265 8250
Phone +1 519 748 6669
Fax +1 519 748 9158
www.bernafon.ca



Amsterdam - Auckland - Berlin - Bern - Brisbane - Copenhagen - Gdańsk - Glasgow - New Jersey - Paris - Rome - Stockholm - Tokyo - Toronto

www.bernafon.com

bernafon[®]
Your hearing • Our passion