




Bertec® Portable Essential

The Bertec® Portable Essential™ provides a variety of assessments, including standard static balance tests, weight shifting, and optional higher-level balance testing for athletes (COBALT™). Static balance and weight shifting activities can be easily integrated into a customized therapy program.

- High-sensitivity dual-balance force plate
- Objective assessment of balance with normative data for comparison
- Convenient for clinical rehabilitation, home health, bedside, and community health programs
- Portable to bring balance services beyond clinic walls
- Training menus for use with biofeedback rehabilitation
- Low profile and light weight

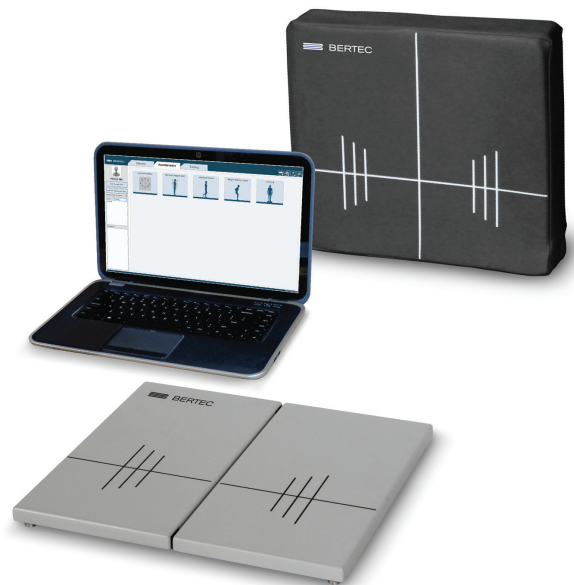

Standard Package

Bertec Balance Advantage Software	Foam balance pad (19" x 20" x 4")
Dual-balance force plate	Database merge and sync
Laptop computer with touchscreen display and remote control	One-year parts and labor with ongoing support
Power cable and USB cable	
Wireless pointing device/mouse	

Assessment Protocols

Modified Clinical Test of Sensory Interaction on Balance (mCTSIB)
Limits of Stability (LOS)
Rhythmic Weight Shift (RWS)
Unilateral Stance (US)
Weight Bearing Squat (WBS)

Featured Training

Quick Training
Mobility
Closed Chain
Weight Shifting
Seated Training

Specifications

High-sensitivity dual-balance force plate
Standardized protocols with normative data
Low 1.5" profile in an 18" x 20" configuration
Modular software design to fit your needs
500 lb load capacity
Patented and HIPAA compliant patient database with merge and sync capabilities across all Bertec® Balance Advantage® products
Portable system powered by USB providing assessments for patients with balance and/or mobility problems

Options

Bertec® Vision Advantage™
COBAL™ (Concussion Balance Test) powered by Bertec®
Carrying case (custom fit to all components)
Extended warranty

For more information, contact Bertec at 614-543-8099 or by email at info@bertec.com