

## ABOUT US

Flex VRA is the first and only visual reinforcement audiometry system that provides the audiologist complete control of the reinforcement stimulus.

Utilizing a wide variety of video sources, Flex VRA is designed to be the most engaging and kid-friendly VRA system available.




"Reliable"  
"Kid-Friendly"  
"Flexible"  
"Simple"  
"Dynamic!"  
"Engaging"



# FLEX

VISUAL REINFORCEMENT AUDIOMETRY

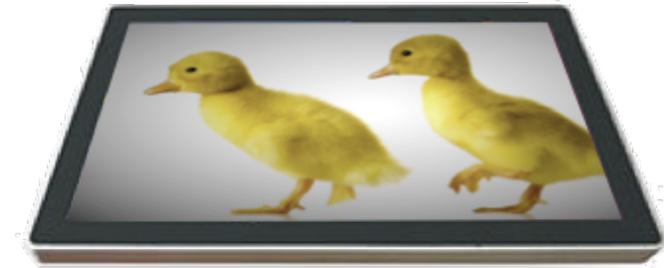

Toy. Animal. Cartoon. One System.

Contact us for information.

### Our Guarantee

#### 90 Day, 100% Money-back Guarantee

Every FLEX VRA purchased from an authorized dealer comes with a 3 year manufacturers warranty, and a 90 day, 100% money-back guarantee. If you're not 100% satisfied with your FLEX VRA system, please let us know and we'll make it right.



(512) 671-0089



100 East Whitestone Boulevard  
Suite 148, #160  
Cedar Park, Texas 78613, United States



sales@flexvra.com

www.flexvra.com

# THE REMOTE

## 1 Stimulus Options

Choose from a wide variety of videos including live-action toy, animal, and cartoon!

## 2 Stimulus Options (with lights)

An animated light bar is available for built-in cartoon, animal, and toy video stimulus.

\*\*\*Animated light bar not available for custom uploaded content.\*\*\*

## 3 ALL and DISTRACT

"ALL" allows for all monitors to be activated simultaneously to localize the patient to the stimulus.

"Distract" will turn the monitors on with one button click and remain on until the button is pressed again. This can assist in distraction during tympanometry or OAE testing.

## 4 Custom Content

Custom content can easily be uploaded to the 3 custom buttons by a provided Flex USB stick.

\*\*\*End User responsible for providing files to upload. Flex VRA will not ship with custom content loaded.\*\*\*

## 5 Volume Control

Volume can be tailored to the patient by the "up" and "down" buttons. In addition, you can also mute the monitors if sound is not desired.

## 6 Video Control

Video skip and repeat options are available for all content, but they are most helpful when navigating custom content.

