


Science **made** smarter

# Caloric Irrigators

Caloric testing  
**made** reliable

Our irrigators offer a choice of stimuli, either water or air, for caloric testing during the VNG test protocol

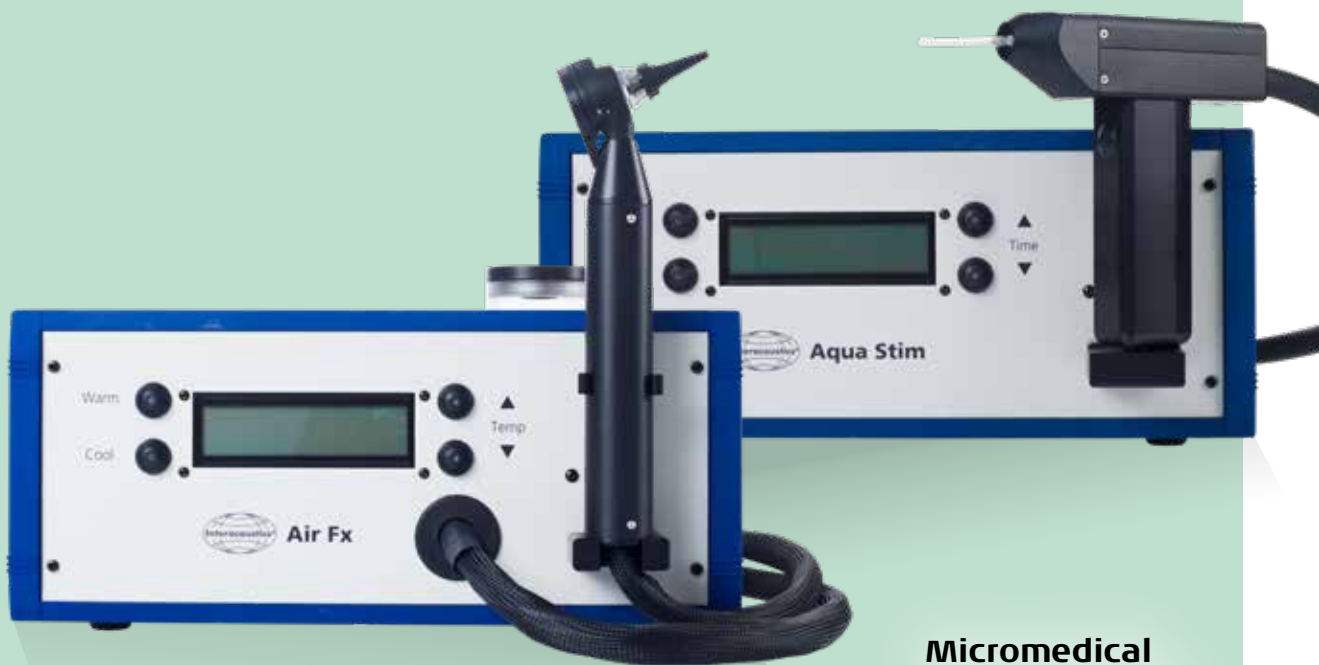

**Micromedical**  
by Interacoustics

  
**Interacoustics**

Audiometry

Tympanometry


ABR

OAE

Hearing Aid Fitting

**Balance**  


# Air Fx Caloric irrigator

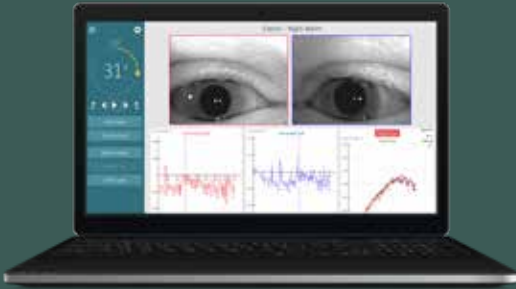

Caloric test with the VisualEyes™515/525 software

Integrated with the  
Interacoustics VNG  
software

## Versatile Caloric Air Irrigator

The AirFx caloric irrigator has many benefits. It offers an otoscopic handle with illumination as well as magnification. In addition, irrigations can be started using the start button on the handle or via the VisualEyes™515/525 or VN415/VO425 software.

The Air Fx is available in 110-130 V and 220-240 V configurations, and is available world-wide.

An integrated water reservoir allows for cool air irrigations down to 20°C regardless of ambient temperature, making it suitable for warm climate use. For warm air irrigations, the Air Fx is capable of heating air up to 50°C.


Temperature Accuracy and Stability is within  $\pm 0,5^{\circ}\text{C}$ , and air flow stability is within  $\pm 0,5$  L/min.

### Highlights:

- Cooling and heating ability: 20°C to 50°C
- Illuminated and magnified otoscope
- Self-contained cooling ability
- Integrated with our VNG software



# Aqua Stim Caloric irrigator

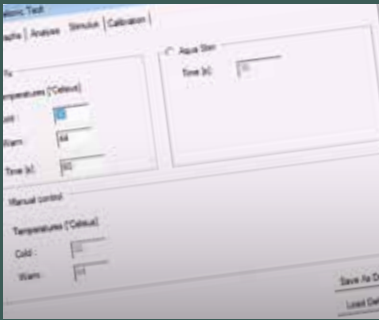

Can be operated from the VisualEyes™515/525 or VN415/VO425 software.

## Versatile Caloric Water Irrigator

The Interacoustics Aqua Stim offers reliable caloric irrigation. The Aqua Stim can be used as a stand-alone device, or you can control it from the VisualEyes™515/525 or VN415/VO425 software.

Irrigation is started by simply pressing the start button on the handle one time to initiate the user-defined caloric protocol.


An external water tank has capacity for approximately 14 irrigations, and eliminates the need for a drain in the examination room.

The Aqua Stim is available world-wide and is supplied as a 110-130 V or 220-240 V model.

### Highlights:

- 30°C cold and 44°C warm irrigations
- Use with our VNG software or as stand-alone
- LED illumination and start button in handle
- Eliminates the need for drain and faucet in examination room

Integrated with the  
Interacoustics VNG  
software

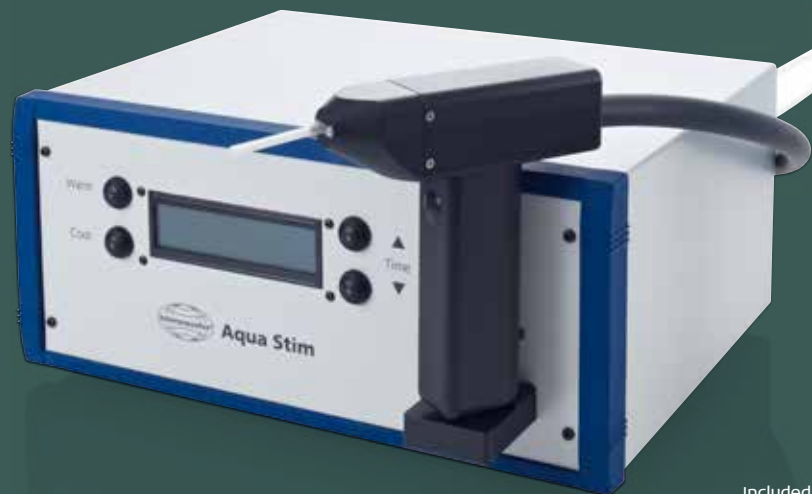

Included external water tank

# Science made smarter

## Interacoustics is more than state-of-the-art solutions

Our mission is clear. We want to lead the way in audiology and balance by translating complexity into clarity:

- Challenges made into clear solutions
- Knowledge made practical
- Invisible medical conditions made tangible and treatable

Our advanced technology and sophisticated solutions ease the lives of healthcare professionals.

We will continue to set the standard for an entire industry. Not for the sake of science. But for the sake of enabling professionals to provide excellent treatment for their millions of patients across the globe.

[Interacoustics-us.com](http://Interacoustics-us.com)

**Interacoustics USA**  
10393 West 70th Street  
Eden Prairie, MN 55344

T +1 800 947 6334  
F +1 952 903 4200


[info@interacoustics-us.com](mailto:info@interacoustics-us.com)  
[interacoustics-us.com](http://interacoustics-us.com)

Go online to  
explore our  
full product  
range

## Related products



**VisualEyes 525**  
Complete VNG solution for  
balance assessment

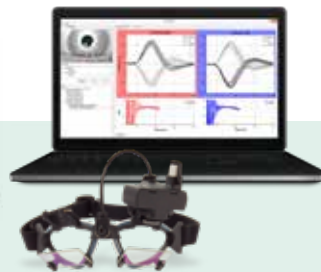

**EyeSeeCam vHIT**  
Video Head Impulse Test



**Orion**  
Rotary Chair for balance  
assessment

### Product specifications

All technical and hardware specifications concerning all products.



**Interacoustics**