

Science made smarter

AA222

Tympanometry
and audiometry
made perfect

The AA222
combines the power
of a clinical middle
ear analyzer with a
complete diagnostic
audiometer.



Interacoustics

Audiometry
~~~~~

Tympanometry  
~~~~~


ABR

OAE

Hearing Aid Fitting

Balance


A unique and powerful diagnostic combo

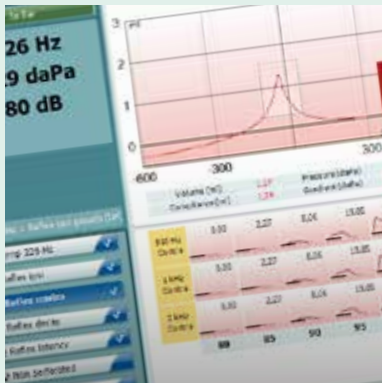

Tympanogram and reflex test results displayed simultaneously



Speech test available in both table and graphic mode

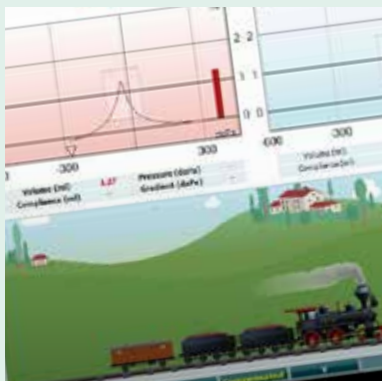

Child distractor - moving choo-choo train



The perfect combination device when space and cost savings need to be considered. The power of a clinical middle ear analyzer with a complete diagnostic audiometer makes the AA222 perfect for a start-up practice or satellite locations.

Comprehensive test selections

- 226 Hz tympanometry (High Hz optional)
- Ipsi/contra reflexes
- Reflex decay
- ETF tests
- Pure tone air & bone conduction
- Speech testing
- Auto HL, Stenger, ABLB, Weber

Outstanding features - impedance

- Child distractor - moving choo-choo train
- Large, 10" color display for easy viewing
- Durable clinical probe
- Optional diagnostic probe for fast screening applications
- Brilliant probe status lights
- Customized test sequences

Outstanding features - audiometry

- Multiple user test setups
- Large, 10" color audiogram display
- Comprehensive library of speech recordings
- Binaural speech testing
- Special QuickSIN test mode
- Operator mic/monitor headset - or -
- Built-in gooseneck mic and speaker

Outstanding performance

- Extremely fast tympanometry and reflex tests
- Customized setups for impedance and audiometry
- Direct printing to external printer
- Included database or NOAH interface
- Diagnostic Suite - merges test data into custom reports

**Tympanometry
& audiometry
in one powerful
device**



Color:	Probe:	Status:
Red		Right ear is selected. Probe is out of the ear.
Blue		Left ear is selected. Probe is out of the ear.
Green		Probe is in the ear and a seal is maintained.
Yellow		Probe is in the ear and is blocked or leaking.
White		The probe has just been attached. Probe status is yet unknown.

Probe status displayed on diagnostic probe / shoulder box for clinical probe



Science made smarter

Interacoustics is more than state-of-the-art solutions

Our mission is clear. We want to lead the way in audiology and balance by translating complexity into clarity:

- Challenges made into clear solutions
- Knowledge made practical
- Invisible medical conditions made tangible and treatable

Our advanced technology and sophisticated solutions ease the lives of healthcare professionals.

We will continue to set the standard for an entire industry. Not for the sake of science. But for the sake of enabling professionals to provide excellent treatment for their millions of patients across the globe.

Interacoustics-us.com

Interacoustics USA

10393 West 70th Street
Eden Prairie, MN 55344

T +1 800 947 6334
F +1 952 903 4200

info@interacoustics-us.com
interacoustics-us.com

Go online to
explore our
full product
range

Related products



AC40
Clinical audiometer



Titan
Middle ear analyzer



Viot™
Video otoscope

Product specifications

All technical and hardware specifications concerning all products can be downloaded from our website.



Interacoustics