



easyScreen

Automated ABR + OAE

Newborn Hearing Screening

Our new combined ABR + OAE Newborn Hearing Screener

With our new easyScreen we unite automated ABR, DPOAE and TEOAE screening capabilities in one device for the first time in MAICO history.

easyScreen is designed to lighten your daily workload in many ways. Small and light-weight, it fits in the palm of your hand. The resistive 4.3" touchscreen makes operation simple even when wearing medical gloves. Achieve accurate results with very fast test times with only a few taps.

easyScreen

Small but powerful - advanced test methods and algorithms achieve accurate results quickly



Use easyScreen's convenient cradle for wireless charging while managing your data.

HearSIM™ - Intuitive software for easy Device and Data Management

Connect the easyScreen to your PC and it's ready for data transfer. Start the transfer with a few simple clicks to store, view and manage patient data in the dedicated HearSIM™ software.

Using HearSIM™ is just as intuitive as operating your easyScreen. Choose to print a full report from HearSIM™ or wirelessly print a test results label directly from easyScreen "on the go" using the optional label printer.



Newborn Hearing Screening has never been so easy

Perfect for hospitals and clinics, easyScreen is fast, easy to handle and space-saving. Even first-time users maneuver easily through the menus and achieve accurate results quickly. The software reminds the screener about important steps before testing and displays information when conditions are not ideal for achieving good quality results.



Powerful automated ABR for fast, accurate results

easyScreen uses the unique CE-Chirp® stimulus with a powerful response detection algorithm. Together they achieve accurate pass or refer results quickly.

The patented CE-Chirp® stimulates all regions of the cochlea at the same time generating a significantly larger ABR than the standard click stimulus used in many other ABR screening systems. This leads to fast, reliable results under normal nursery conditions.

Test both ears simultaneously to save valuable time



easyScreen Cart


Our optional easyScreen rolling cart keeps the components secure and conveniently stores all necessary supplies for easy transport from room to room.

The greatest care for the tiniest ears

easyScreen's OAE probe is lightweight with comfortable eartip options for the best fit in infant ear canals. Our ABR supports use of eartips or EarCups with our insert earphones or you can use the OAE probe for ABR too. Choose the option that suits your program preferences and budget.



Sanibel
Supply

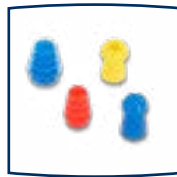

Standard Components



easyScreen device



Cradle



Preamplifier



Insert earphones



OAE Probe



Sanibel eartip starter kit



HearSIM™ Software

Optional Accessories



Wireless label printer



Cart



MAICO Diagnostics

10393 West 70th Street
Eden Prairie, MN 55344

Tel.: 888.941.4201
Fax: 952.903.4100

e-mail: info@maico-diagnostics.com
Web: www.maico-diagnostics.com