

NEW



## MedRx WREM

Wireless REM/LSM System

## Wireless REM


We've Expanded our  
Wireless Portfolio to Include  
a REM Only Portion.

REM AutoFit Compatible.

# MedRx<sup>®</sup>

## 888-392-1234

To Schedule An Online Demo

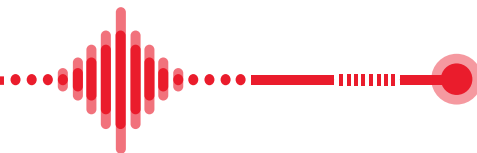




## WREM Standard Accessories

- Wireless Probe Mics with Bow and Charging Base
- Monitor Headphone
- Speaker
- Patient Response Switch
- USB-C Cable
- External Power Supply
- Probe Tube Pack
- Software & Manuals
- Carrying Case

## REM Optional Accessories

- RECD Coupler


## REM

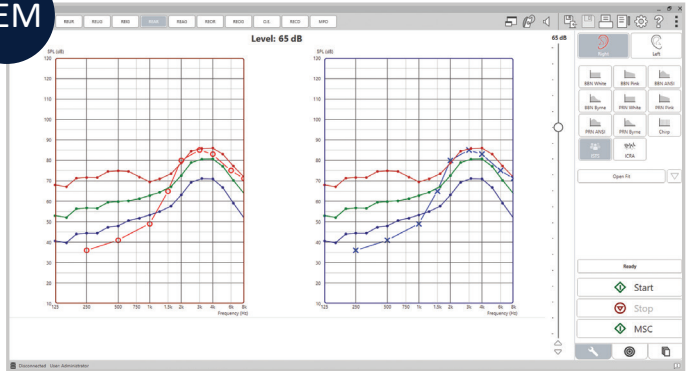

### Real Ear Measurements

IEC/ANSI standard REM functions are provided for: REUR, REUG, REIG, REAG, REAR, REOR, REOG and RECD.


## LSM

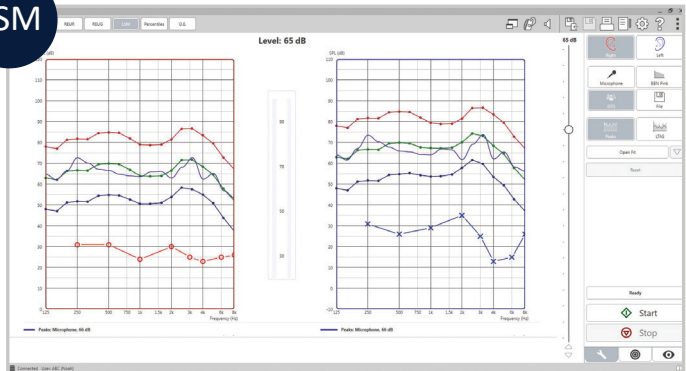

### Live Speech Mapping

Live Speech Mapping (LSM) a quick and successful approach to hearing instrument fitting.

## HLS MHA

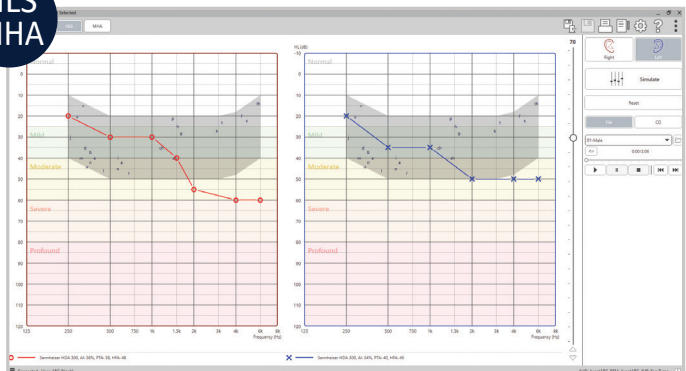

### HLS/MHA

Easily simulate hearing loss and hearing instruments for patients.



Backside of WREM  
Ports for: Monitor, Speaker (Line Out),  
Speakers - Right & Left, USB-C and Power

# MedRx Wireless REM/LSM System

- Binaural live speech mapping and real ear measurements with wireless probe mics
- High-frequency wireless fittings up to 12,500 Hz
- Studio software with Probe Tube Fitting Guide
- REM Autofit Capable
- The REM Software includes targets for MSS (Modified Speech Spectrum), DSL v5.0, NAL-NL1 and NAL-NL2
- 3rd Party Counseling and Demonstration with Hearing Loss Simulator and Master Hearing Aid
- PC-Powered with USB-C Connection and Set-up
- Portable and Compact at Only 20 cm x 6 cm x 5 cm (L x W x H)
- Noah, TIMS, Blueprint OMS, Sycle and OtoAccess 2 Compatible



## MedRx WREM

MedRx is excited to launch the next product in our portfolio of wireless devices, the Wireless REM. Perhaps the Wireless REM's most notable feature is its new wireless probe microphones. This enhancement will minimize fumbling around unruly wires and lead to greater efficiency and maneuverability in your sessions. The wireless probe microphones automatically connect to your device through a private connection used only to communicate between the two instruments, limiting interference and connection barriers.

## Wireless Real Ear Measurement and Live Speech Mapping

MedRx's modern REM and LSM software modules allow clinicians to fit hearing aids more accurately than ever, improving patient satisfaction and reducing hearing aid returns. Powered only by the USB on your computer, the WREM allows clinicians to fit up to 12,500 Hz. The WREM also uses two wireless Bluetooth probe microphones that connect to your system automatically for each fitting session.

To reach 12,500Hz, we use high end wireless transmissions to allow for added bandwidth making our probes high-frequency capability. This makes the Wireless REM suitable for patients who need a unique gain on their devices or are using hearing aids with more advanced high-frequency outputs.

As with all MedRx products, the Wireless REM is 100% computer-based, ideal for a desktop or a laptop. Our real ear and live speech mapping modules within our Studio software are included with each Wireless REM unit and were intentionally engineered for simple, straightforward navigation and functionality.

The included software also features a host of handy tools (including the MHA and HLS) and customizable settings and views. One tool sure to expedite fittings is our Probe Tube Depth Guide, this new guide allows clinicians to insert probe microphone tubes faster and with more confidence during fittings—taking the fear out of probe tube placements for clinicians and patients.

# MedRx WREM

## Technical Specs

**Standards:** IEC 60601-1:2 Class II, IEC 60601-1-2 Class A, IEC 60645-1, Medical Device Regulation

**REM Speech+ Standards:** ANSI S3.46, IEC 61669, EN 61669

**Dimensions:** Approx. 18 cm x 6 cm x 5 cm (L x W x H)  
Approx. 7" x 2" x 1.5" (L x W x H).

**Weight:** < 500 g < 1 lb

**System Modality:** Real Ear Measurement; Binaural Live Speech Mapping; Hearing Loss Simulator; Hearing Instrument Simulator.

**Probe Microphones (L/R):** Dual Electret Microphone Elements (2 Probe Microphones)

**Probe Microphone Tube:** Silicone 1.0 mm Nominal Diameter

**Measurement Frequency Range:** 125 – 12500 Hz

**Test Stimuli:** Broadband Noise and Synthesized Random Noise – Pink, White, Byrne LTASS and ANSI weighted; ICRA; ISTS; Microphone, File, CD-ROM for Live Speech Mapping, Chirp

**Test Stimulus Levels at 1m:** 45 – 90 dB SPL in 1 dB Steps – 200 Hz thru 12500 Hz (depending on speaker wattage & efficiency)

**Test Stimulus Accuracy:** ± 3 dB SPL

**Analysis Mode:** User Selectable 1/3, 1/6, 1/12, 1/24, 1/48 Octave Bands

**ANSI S3.46 Test Available IEC/EN 61669:** Real Ear Unaided Response, Real Ear Unaided Gain; Real Ear Insertion Gain; Real Ear Occluded Response; Real Ear Occluded Gain; Real Ear Aided Response; Real Ear Aided Gain

**Other Test Available:** Live Speech Mapping with Peaks and LTASS analysis; Real Ear to Coupler Difference, Occlusion Effect, Percentile Analysis

**Prescription Methods:** NAL-RP; 1/3 Gain; 1/2 Gain; Berger; Pogo 1; Pogo 2; FIG6; DSL m[I/O] NAL-NL1; NAL-NL2

**External Connections:** Power Connection USB 3.0 Input 5.0 Volt Bus; Line-Output Jack (Speakers) 3.5 mm Stereo Jack; Speaker Output (Internal Amplifier) (2) 3.81mm Pluggable Spring Clamp; Monitor Headset Jack 3.5 mm Stereo Jack; Power Jack 2.1 mm X 5.5 mm.

**Data Connection:** USB

**Mode of Operation:** Continuous

**Warm up Time:** Less than 5 min after USB connection

**Power Connection:** USB 3.0 Input 5.0 Volt Bus

**USB 3.0 Input:** Standard USB "C" Socket

**Power Consumption:** Less than 1800 mA at 15 VDC / less than 900mA at 5 VDC

**Operating Temperature:** 10°C to 35°C

**Operating Humidity:** 30% to 90%

**Storage Temperature:** -20°C to 50°C

**Storage Humidity:** 10% to 90%

### MedRx Minimum Computer Specs:

Windows® PC compatible computer, Intel™ i5, 2.0 GHz or better. 4 GB RAM. 20 GB free hard drive space. Available 2.0 USB Port. Windows 10 or 11 Professional (32 Or 64-Bit).

# MedRx®

Good Things Come In Small Packages

### MedRx, Inc

1200 Starkey Rd. #105

Largo, FL 33771 U.S.A.

888-392-1234 • 727-584-9600

Fax: 727-584-9602

Email: [sales@medrx-diagnostics.com](mailto:sales@medrx-diagnostics.com)

[www.medrx-diagnostics.com](http://www.medrx-diagnostics.com)

MX-WREM-DS-1 Rev. 1

The probe microphone's rechargeable battery provides 6 hours of continuous use

