

OTICON | Zircon

Technical data sheet

miniRITE T

60 85 100 105



	Zircon 1	Zircon 2	
Speech Understanding	OpenSound Navigator™	•	-
	- Balancing power effect	40%	-
	- Max. noise removal difficult/simple	6 dB / 0 dB	-
	Multiband Adaptive Directionality	-	•
	Noise Reduction	-	•
	Speech Guard™	•	-
	Single Compression	-	•
Sound Quality	Frequency lowering	Speech Rescue™	Speech Rescue™
	Fitting Bandwidth*	8 kHz	8 kHz
	Bass Boost (streaming)	•	•
Listening Comfort	Processing Channels	48	48
	Feedback Management	SuperShield & Feedback shield	SuperShield & Feedback shield
	Transient Noise Management	On/Off	-
Personalization & Optimizing Fitting	Wind Noise Management	•	•
	Fitting Bands	14	12
	Multiple Directionality options	•	•
	Adaptation Management	•	•
	Oticon Firmware Updater	•	•
Connecting to the world	Fitting Formulas	NAL-NL1/NAL-NL2, DSL 5.0	NAL-NL1/NAL-NL2, DSL 5.0
	Hands-free communication**	•	•
	Direct streaming***	•	•
	Oticon ON app & Oticon RemoteCare app	•	•
	ConnectClip	•	•
	EduMic	•	•
	Remote Control 3.0	•	•
	TV Adapter 3.0	•	•
	Phone Adapter 2.0	•	•
	Tinnitus SoundSupport™	•	•
CROS/BiCROS support	•	•	

*Bandwidth accessible for gain adjustments during fitting

**Available for Oticon Zircon from FW 1.1 with select iPhone models

***From iPhone®, iPad®, iPod touch®, and select Android™ devices

Operating Conditions

Temperature: +1°C to +40°C (34°F to 104°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Storage and transportation conditions

Temperature and humidity should not exceed the below limits for extended periods during transportation and storage.

Transportation

Temperature: -25°C to +60°C (-13°F to 140°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Storage

Temperature: -25°C to +60°C (-13°F to 140°F)

Humidity: 5% to 93% relative humidity, non-condensing

Atmospheric pressure: 700 hPa to 1060 hPa

Apple, the Apple logo, iPhone, iPad, and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

Oticon Zircon miniRITE T offers a discreet design with LED-light to make handling easy. The style features telecoil and a double push-button. It is a Made for iPhone® hearing aid and compatible with the new Android protocol for Audio Streaming for Hearing Aids (ASHA) - making it possible to stream directly from iPhone, iPad®, iPod touch® and select Android™ devices.

OpenSound Navigator™ provides access to speech in 360° making the listener more easily aware of what is going on in the surroundings.

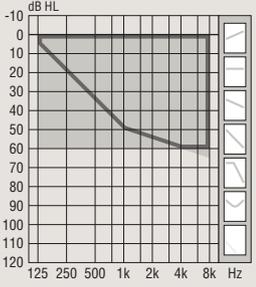
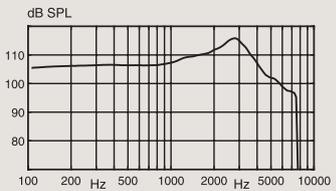
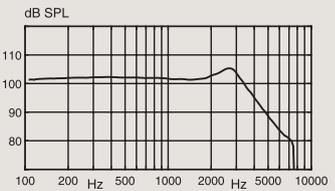
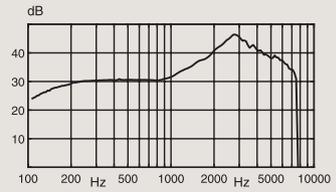
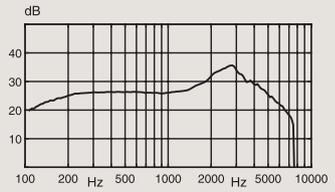
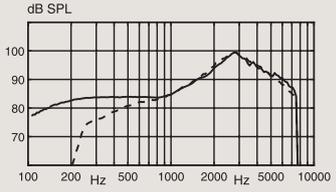
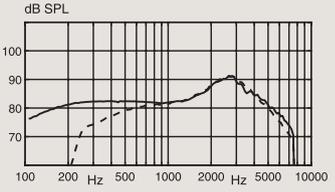
Speech Guard™ provides more natural and clear speech sounds making the details in speech stand out more.

The Polaris™ platform provides a tremendous speed and memory capacity for audiological processing and connectivity options. New features can be added and updates performed wirelessly.

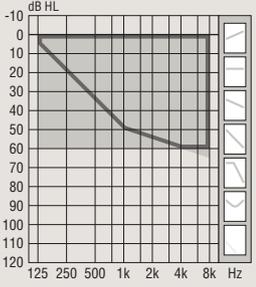
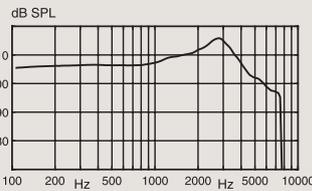
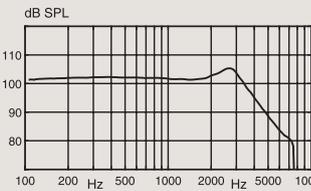
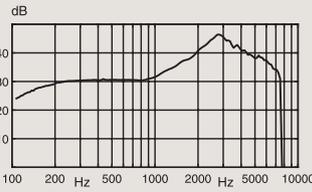
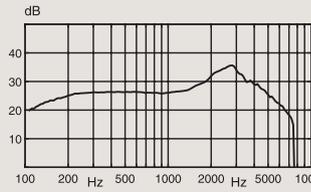
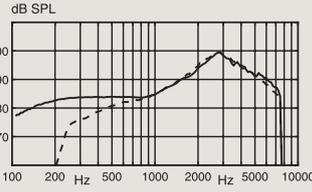
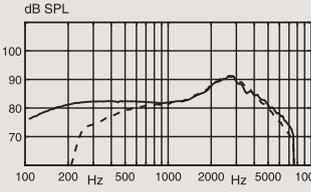


For information on compatibility, please visit www.oticon.com/support/compatibility

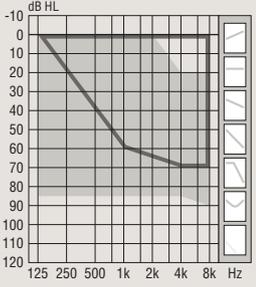
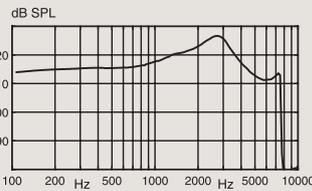
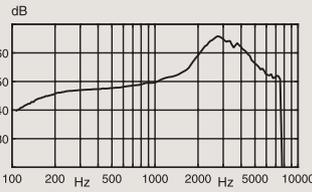
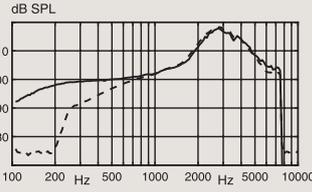
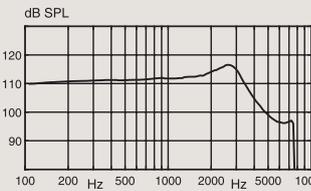
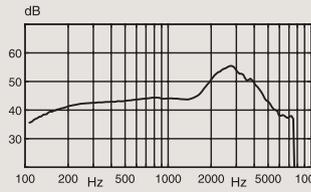
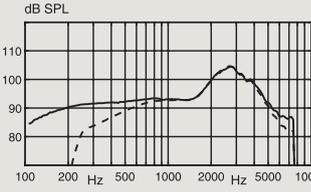
oticon
life-changing technology

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 60		OSPL90 	OSPL90 
		Full-on gain 	Full-on gain 
Technical information Omnidirectional mode is used unless otherwise stated.		Frequency response 	Frequency response 
— Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m			
OSPL90	Peak 1600 Hz HFA-OSPL90	116 dB SPL 110 dB SPL 111 dB SPL	105 dB SPL 102 dB SPL 103 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	46 dB 37 dB 38 dB	36 dB 29 dB 30 dB
Reference test gain		30 dB	26 dB
Frequency range		100-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	68 dB SPL 88 dB SPL -	- - 85/85 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 2 % < 3 % < 2 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	18 dB SPL 26 dB SPL	16 dB SPL 27 dB SPL
Battery consumption ²	Typical Quiescent	2.2 mA 2.2 mA	2.2 mA 2.2 mA
Battery life, artificial measurement, hours ³		80	80
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		55-60	

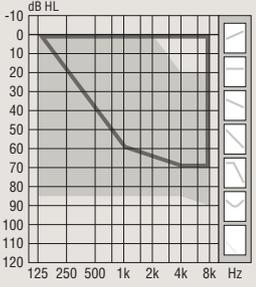
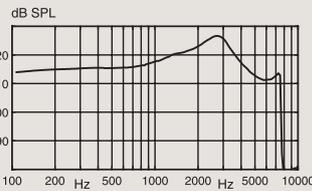
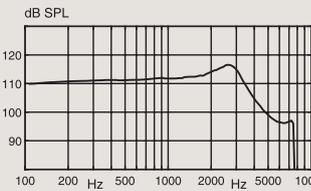
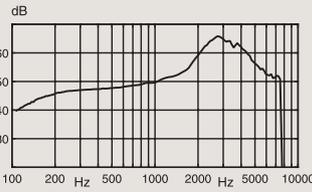
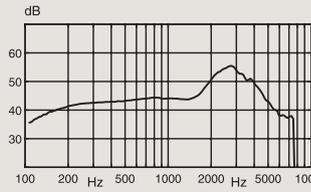
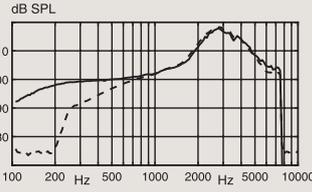
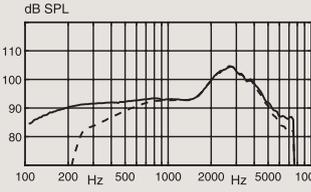
1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 <p>60</p> <p>Legend: <input checked="" type="checkbox"/> Mold, Bass & Power dome <input type="checkbox"/> OpenBass dome</p> <p>Technical information Omnidirectional mode is used unless otherwise stated.</p>		<p>OSPL90</p> 	<p>OSPL90</p> 
		<p>Full-on gain</p> 	<p>Full-on gain</p> 
		<p>Frequency response</p>  <p>— Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m</p>	<p>Frequency response</p> 
OSPL90	Peak 1600 Hz HFA-OSPL90	116 dB SPL 110 dB SPL 111 dB SPL	105 dB SPL 102 dB SPL 103 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	46 dB 37 dB 38 dB	36 dB 29 dB 30 dB
Reference test gain		30 dB	26 dB
Frequency range		100-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	68 dB SPL 88 dB SPL -	- - 85/85 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 2 % < 3 % < 2 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	18 dB SPL 26 dB SPL	16 dB SPL 27 dB SPL
Battery consumption ²	Typical Quiescent	2.2 mA 2.2 mA	2.2 mA 2.2 mA
Battery life, artificial measurement, hours ³		80	80
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		55-60	

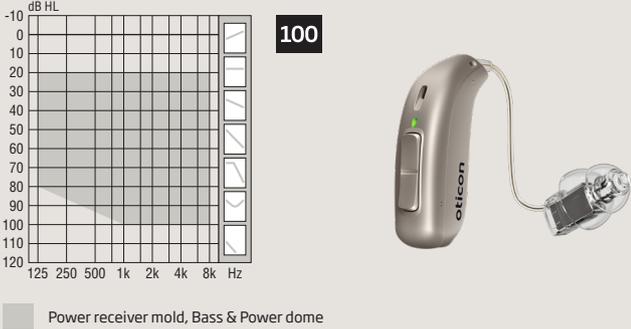
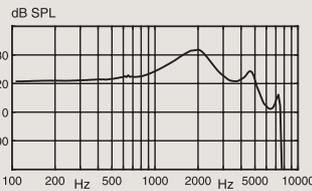
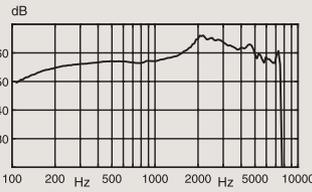
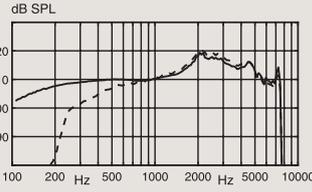
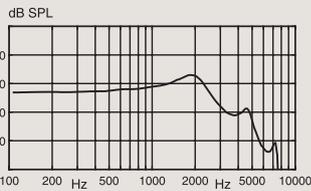
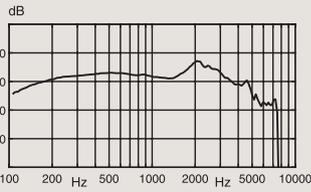
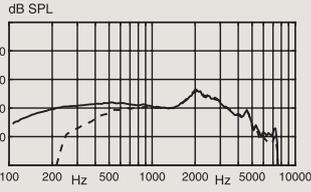
1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006	
 <p>85</p>  <p>Legend: <input checked="" type="checkbox"/> Mold, Bass & Power dome <input type="checkbox"/> OpenBass dome</p> <p>Technical information Omnidirectional mode is used unless otherwise stated.</p>		<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p>  <p>— Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m</p>	<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p> 	
	OSPL90	Peak 1600 Hz HFA-OSPL90	127 dB SPL 121 dB SPL 122 dB SPL	117 dB SPL 113 dB SPL 114 dB SPL
	Full-on gain ¹	Peak 1600 Hz HFA-FOG	66 dB 53 dB 56 dB	55 dB 45 dB 48 dB
	Reference test gain		46 dB	37 dB
Frequency range		100-7500 Hz	100-7500 Hz	
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	84 dB SPL 104 dB SPL -	- - 96/96 dB SPL	
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 2 % < 4 % < 5 %	< 2 % < 2 % < 2 %	
Equivalent input noise level	Omni Dir	21 dB SPL 28 dB SPL	17 dB SPL 27 dB SPL	
Battery consumption ²	Typical Quiescent	2.3 mA 2.2 mA	2.4 mA 2.2 mA	
Battery life, artificial measurement, hours ³		75	75	
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-60		

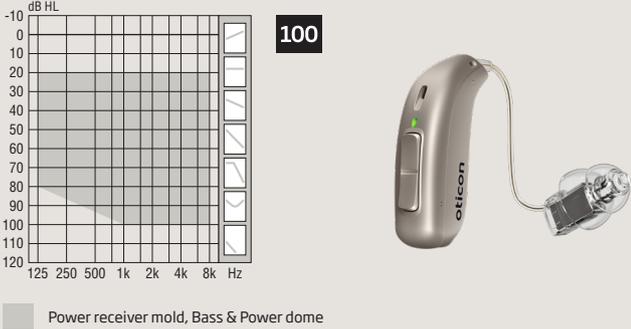
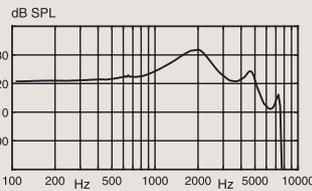
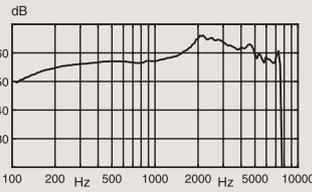
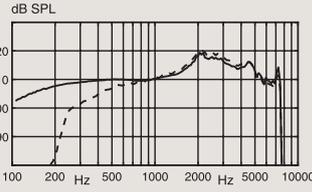
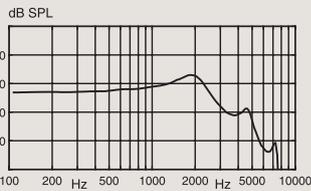
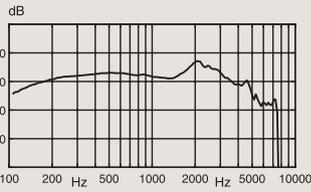
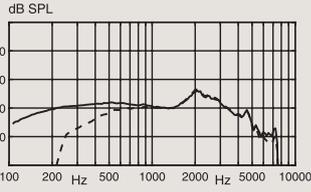
1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 Mold, Bass & Power dome OpenBass dome Technical information Omnidirectional mode is used unless otherwise stated.		OSPL90 	OSPL90 
		Full-on gain 	Full-on gain 
		Frequency response 	Frequency response 
		— Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m	
OSPL90	Peak 1600 Hz HFA-OSPL90	127 dB SPL 121 dB SPL 122 dB SPL	117 dB SPL 113 dB SPL 114 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	66 dB 53 dB 56 dB	55 dB 45 dB 48 dB
Reference test gain		46 dB	37 dB
Frequency range		100-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	84 dB SPL 104 dB SPL -	- - 96/96 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	< 2 % < 4 % < 5 %	< 2 % < 2 % < 2 %
Equivalent input noise level	Omni Dir	21 dB SPL 28 dB SPL	17 dB SPL 27 dB SPL
Battery consumption ²	Typical Quiescent	2.3 mA 2.2 mA	2.4 mA 2.2 mA
Battery life, artificial measurement, hours ³		75	75
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-60	

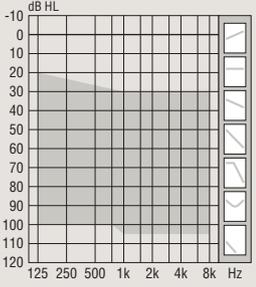
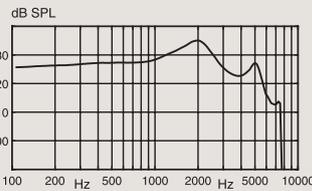
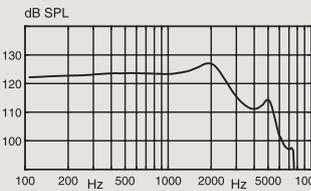
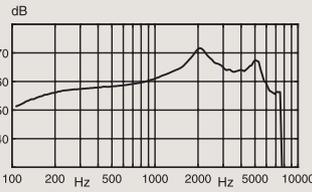
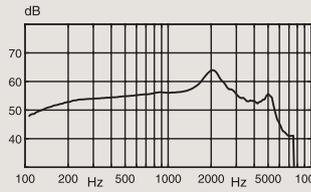
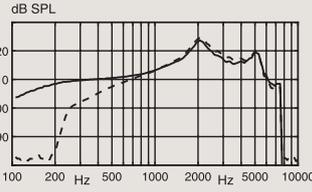
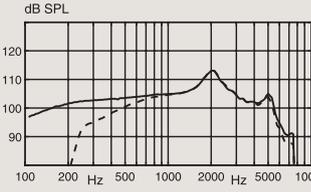
1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006	
 <p>100</p> <p>Power receiver mold, Bass & Power dome</p> <p>Technical information Omnidirectional mode is used unless otherwise stated.</p> <p>Warning to the hearing aid dispenser The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the hearing aid, as there may be risk of impairing the remaining hearing of the hearing aid user.</p> <p>— Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m</p>		<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p> 	<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p> 	
	OSPL90	Peak 1600 Hz HFA-OSPL90	132 dB SPL 130 dB SPL 127 dB SPL	123 dB SPL 122 dB SPL 119 dB SPL
	Full-on gain ¹	Peak 1600 Hz HFA-FOG	66 dB 60 dB 61 dB	57 dB 53 dB 53 dB
	Reference test gain		53 dB	42 dB
Frequency range		100-7500 Hz	100-7500 Hz	
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	91 dB SPL 111 dB SPL -	- - 101/101 dB SPL	
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	<9 % <6 % <3 %	<2 % <2 % <2 %	
Equivalent input noise level	Omni Dir	16 dB SPL 25 dB SPL	16 dB SPL 28 dB SPL	
Battery consumption ²	Typical Quiescent	2.2 mA 2.2 mA	2.3 mA 2.2 mA	
Battery life, artificial measurement, hours ³		80	75	
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-60		

1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 <p>100</p> <p>Power receiver mold, Bass & Power dome</p>		<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p> 	<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p> 
	<p>— Acoustic input: 60 dB SPL</p> <p>- - - Magnetic input: 31.6 mA/m</p>		
OSPL90	Peak 1600 Hz HFA-OSPL90	132 dB SPL 130 dB SPL 127 dB SPL	123 dB SPL 122 dB SPL 119 dB SPL
Full-on gain ¹	Peak 1600 Hz HFA-FOG	66 dB 60 dB 61 dB	57 dB 53 dB 53 dB
Reference test gain		53 dB	42 dB
Frequency range		100-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field 10 mA/m field SPLITS L/R	91 dB SPL 111 dB SPL -	- - 101/101 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz 800 Hz 1600 Hz	<9 % <6 % <3 %	<2 % <2 % <2 %
Equivalent input noise level	Omni Dir	16 dB SPL 25 dB SPL	16 dB SPL 28 dB SPL
Battery consumption ²	Typical Quiescent	2.2 mA 2.2 mA	2.3 mA 2.2 mA
Battery life, artificial measurement, hours ³		80	75
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-60	

1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

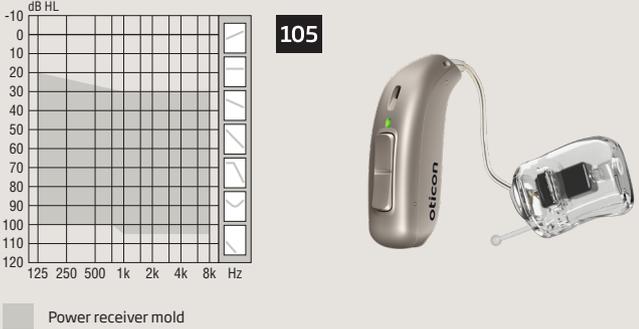
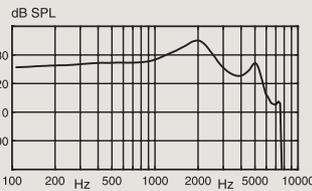
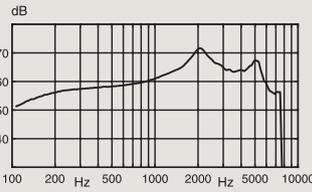
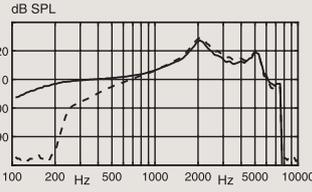
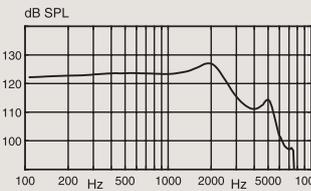
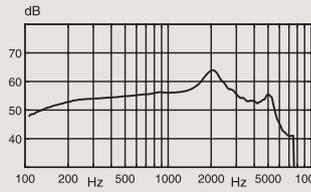
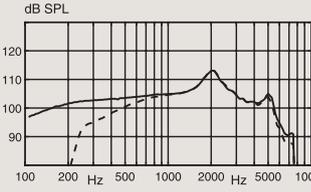
		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 <p>Power receiver mold</p> <p>Technical information Omnidirectional mode is used unless otherwise stated.</p> <p>Warning to the hearing aid dispenser The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the hearing aid, as there may be risk of impairing the remaining hearing of the hearing aid user.</p> <p>— Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m</p>		OSPL90 	OSPL90 
		Full-on gain 	Full-on gain 
		Frequency response 	Frequency response 
OSPL90	Peak	135 dB SPL	127 dB SPL
	1600 Hz	133 dB SPL	126 dB SPL
	HFA-OSPL90	131 dB SPL	123 dB SPL
Full-on gain ¹	Peak	72 dB	64 dB
	1600 Hz	66 dB	59 dB
	HFA-FOG	65 dB	58 dB
Reference test gain		58 dB	47 dB
Frequency range		100-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field	96 dB SPL	-
	10 mA/m field	116 dB SPL	-
	SPLITS L/R	-	106/106 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz	< 4 %	< 2 %
	800 Hz	< 4 %	< 2 %
	1600 Hz	< 4 %	< 2 %
Equivalent input noise level	Omni	15 dB SPL	16 dB SPL
	Dir	24 dB SPL	27 dB SPL
Battery consumption ²	Typical	2.3 mA	2.4 mA
	Quiescent	2.2 mA	2.2 mA
Battery life, artificial measurement, hours ³		80	75
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-60	

1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.

2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.

3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.

4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

		Ear Simulator Measured according to IEC 60118-0:1983/AMD1:1994, IEC 60118-0:2015, IEC 60118-1:1995+AMD1:1998 CSV and IEC 60318-4:2010	2CC Coupler Measured according to ANSI S3.22-2014, IEC 60118-0:2015 and IEC 60318-5:2006
 <p>105</p> <p>Power receiver mold</p> <p>Technical information Omnidirectional mode is used unless otherwise stated.</p> <p>Warning to the hearing aid dispenser The maximum output capability of the hearing aid may exceed 132 dB SPL (IEC 711). Special care should be exercised in selecting and fitting the hearing aid, as there may be risk of impairing the remaining hearing of the hearing aid user.</p> <p>— Acoustic input: 60 dB SPL - - - Magnetic input: 31.6 mA/m</p>		<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p> 	<p>OSPL90</p>  <p>Full-on gain</p>  <p>Frequency response</p> 
OSPL90	Peak	135 dB SPL	127 dB SPL
	1600 Hz	133 dB SPL	126 dB SPL
	HFA-OSPL90	131 dB SPL	123 dB SPL
Full-on gain ¹	Peak	72 dB	64 dB
	1600 Hz	66 dB	59 dB
	HFA-FOG	65 dB	58 dB
Reference test gain		58 dB	47 dB
Frequency range		100-7500 Hz	100-7500 Hz
Telecoil output (1600 Hz)	1 mA/m field	96 dB SPL	-
	10 mA/m field	116 dB SPL	-
	SPLITS L/R	-	106/106 dB SPL
Total harmonic distortion (Input 70 dB SPL)	500 Hz	< 4 %	< 2 %
	800 Hz	< 4 %	< 2 %
	1600 Hz	< 4 %	< 2 %
Equivalent input noise level	Omni	15 dB SPL	16 dB SPL
	Dir	24 dB SPL	27 dB SPL
Battery consumption ²	Typical	2.3 mA	2.4 mA
	Quiescent	2.2 mA	2.2 mA
Battery life, artificial measurement, hours ³		80	75
Expected battery life, hours (battery size 312 - IEC PR41) ⁴		50-60	

1) Measured with the gain control of the hearing aids set to their full-on position minus 20 dB and with an input SPL of 70 dB. This is to obtain a gain response equal to the full-on gain response from e.g. IEC 60118-0:1983+A1:1994 but without influence of feedback.
 2) Battery current is measured according to IEC 60118-0:1983/AMD1:1994 §7.11, IEC 60118-0:2015 §7.7 and ANSI S3.22:2014 §6.13 after a settling time of minimum 3 minutes.
 3) Based on the standardized battery consumption measurement (IEC 60118-0:1983/AMD1:1994). The actual battery life depends on battery quality, use pattern, active feature set, hearing loss and sound environment.
 4) Real usage battery life is shown as an estimated interval based on mixed use cases with variable amplification settings and variable input levels, incl. direct stereo streaming from a TV (25% of the time) and streaming from a mobile phone (6% of the time).

Headquarters
Oticon A/S
Kongebakken 9
DK-2765 Smørum
Denmark



SBO Hearing A/S
Kongebakken 9
DK-2765 Smørum
Denmark